

TUNICS-PLUS

USER MANUAL



Tunable External Cavity Laser

Legal Notice

TUNICS-Plus Tunable External Cavity Laser Instruction Manual

This documentation is provided as an instruction manual to Yenista customers and potential customers only.

Copyright Notice

Copyright © 1999-2002 by Yenista.

The contents of this publication may not be reproduced in any part or as a whole, transcribed, stored in a retrieval system, translated into any language, or transmitted in any form or by any means, electronic, mechanical, magnetic, optical, chemical, photocopying, manual, or otherwise, without the prior written permission of Yenista.

Published by Yenista. All rights reserved.

Trademarks

Yenista, TUNICS, and Yenista products are registered trademarks of Yenista. Other trademarks mentioned in this publication are used for identification purposes only and may be properties of their respective owners.

Yenista US Patents

The proprietary configuration of TUNICS is patented under US patents # 5,594,744, # 5,802,085 and # 6,252,718 B1.

Disclaimers

No Warranties: This documentation is provided “as is” without any express or implied warranty of any kind.

Limitation of Liability: Yenista does not assume any liability arising out of the application or use of any products, or software described herein. In no event shall Yenista or its affiliate companies be liable for any damages whatsoever (including, without limitation, consequential or incidental damages, damages for loss of profits, or otherwise) arising out of the use of the information provided in this documentation.

The contents of this publication have been carefully checked for accuracy. However, Yenista makes no warranty as to the completeness, correctness, or accuracy of the information contained within this documentation. In the interest of continued product improvement, Yenista further reserves the right to make changes in any products described herein without notice. This publication is subject to change without notice.

Contents

LEGAL NOTICE.....	I
CONTENTS.....	I
LIST OF TABLES	V
LIST OF FIGURES	VII
 ABOUT THIS MANUAL	 3
PURPOSE.....	3
AUDIENCE	3
PREREQUISITES	3
ABOUT THIS MANUAL	4
INTRODUCTION	4
CHAPTERS/APPENDICES.....	4
GENERAL SAFETY CONSIDERATIONS.....	5
<i>Terms and Symbols in this Manual</i>	<i>5</i>
<i>Terms on the Product</i>	<i>5</i>
<i>Symbols and Labels on the Product</i>	<i>5</i>
<i>Laser Safety.....</i>	<i>6</i>
General Laser Specifications	6
Laser Class.....	6
Absolute Maximum Ratings	6
<i>Equipment Safety Warning.....</i>	<i>7</i>
<i>Electrical Safety</i>	<i>7</i>
<i>Ventilation of the Unit.....</i>	<i>7</i>
CONTACTING YENISTA	8
 GENERAL INFORMATION	 11
PRODUCT OVERVIEW	11
OPTICAL DESIGN.....	12
BENEFITS.....	13
SPECIFICATIONS	15
 GETTING STARTED	 19
UNPACKING.....	19
STANDARD ACCESSORIES	19
OPTIONS AND ORDERING INFORMATION	19
<i>Models</i>	<i>19</i>
<i>Options</i>	<i>20</i>
<i>Accessories.....</i>	<i>20</i>
<i>Ordering code.....</i>	<i>20</i>
INSTALLATION & FIRST-TIME OPERATION	21
<i>Installation.....</i>	<i>21</i>

Contents

<i>Initializing the Laser</i>	22
<i>Enabling the Laser Output</i>	22
<i>Optical Connection to the Unit</i>	22
<i>Shutting Down the TUNICS-Plus Laser</i>	23
FRONT & REAR PANELS	27
FRONT PANEL.....	27
<i>Description of the Front Panel Functions Keys</i>	27
<i>DATA Area</i>	27
<i>MODE Area</i>	28
<i>SETUP Area</i>	29
<i>General Parameter Setting Principle</i>	29
Using the Numeric Keypad	29
Using the Rotating Modify Knob	30
REAR PANEL	31
<i>Description of the Rear Panel Connectors</i>	31
Monitor (Scaled) Outputs.....	31
TTL Outputs.....	31
Inputs.....	31
Remote Interfaces.....	31
OPERATION	35
INTRODUCTION	35
MOST COMMONLY USED FUNCTIONS.....	36
<i>To Change the Emission Wavelength</i>	36
<i>To Change the Wavelength Step</i>	36
<i>To Change the Optical Output Power</i>	36
<i>To Change the Diode Current</i>	37
<i>To Switch Between Constant-Power and Constant-Current Modes</i>	37
<i>To Use the Fine Scanning Mode</i>	37
ADVANCED FUNCTIONS	38
<i>To Implement Step-by-Step Wavelength Scanning</i>	38
<i>To Implement Continuous Wavelength Sweep (hop-free)</i>	39
Wavelength sweep commands on front-panel	39
Sweep speed	40
<i>To Perform Power Calibration</i>	40
<i>To Turn On Backlash Suppression of the Micrometer Screw</i>	42
<i>To Change the GPIB Address</i>	42
<i>To Turn On/Off the RS-232 Remote Control Mode</i>	42
<i>To Turn On Coherence Control</i>	42
INTERNAL WAVELENGTH REFERENCING.....	43
<i>Self-calibration at start-up</i>	43
<i>User-command for auto-calibration</i>	43
ACTIVE CAVITY CONTROL.....	45
<i>Active cavity control operating requirements</i>	45
<i>To disable/enable active cavity control</i>	45
ACTIVE WAVELENGTH MONITORING	46
<i>Instrument Setup and System Requirements</i>	46
<i>Wavelength Monitoring Scheme</i>	48
<i>How to Perform Wavelength Monitoring</i>	48
AUXILIARY INPUTS AND OUTPUTS	50
<i>External Modulation</i>	50

Intensity Modulation - Low Frequency	50
Low Frequency Modulation Setup	50
Using Low Frequency Modulation	51
Resulting LF-Modulated Diode Current	51
Intensity Modulation - High Frequency	52
High Frequency Modulation Setup	52
Using High Frequency Modulation	53
Resulting HF-Modulated Diode Current	54
<i>Internal Modulation</i>	55
<i>Synchronization Signal</i>	55
<i>Wavelength Modulation</i>	55
Wavelength fine scanning	55
Coherence control	55
<i>Monitor Outputs</i>	55
<i>Remote Interlock</i>	56
REMOTE CONTROL	59
INTRODUCTION	59
RS-232 C REMOTE CONTROL	59
Physical Interface	59
Message Format	60
Response Messages	61
Format of Numerics	61
Changing Parameter Settings	62
Changing the Operating Mode	63
Reading System Parameters	64
Other System Commands	65
IEEE-488.1 REMOTE CONTROL	66
GPIB Overview	66
GBIP Functions Available on TUNICS-Plus	66
Checking System Status with the Status Word Variable	67
Definition of the Status Word	67
Use of the Service Request (SRQ)	68
Syntax for Parameter Commands	69
Syntax and Message Conventions	69
Response messages	69
Value format	69
Synchronization of Messages	70
Error Handling	71
Changing the GPIB Address	71
Local Lockout	72
IEEE-488.1 Control Commands	72
IEEE Standard Commands	72
Changing Parameter Settings	72
Changing Operating Mode	73
System and Parameter Queries	74
SERVICE INFORMATION	77
PERFORMANCE VERIFICATION	77
Equipment Required	77
Preparation	77
Initializing the System	78
Connecting the Laser Output	78
Enabling the Laser Output	78
Shutting Down the Laser	79

Contents

<i>Absolute Wavelength Accuracy</i>	79
<i>Wavelength Range and Output Power</i>	79
<i>Side Mode Suppression</i>	80
<i>Wavelength and Power Calibration</i>	80
MAINTENANCE	81
<i>Cleaning and Caring for Optical Connectors and Fibers</i>	81
<i>Cleaning of the Instrument</i>	82
<i>User Calibration</i>	82
<i>Replacing the Line Fuse</i>	82
<i>Packaging for Shipment</i>	82
APPENDICES	85
APPENDIX A: ERROR AND WARNING MESSAGES	85
APPENDIX B: VALID RANGE FOR PARAMETER SETTINGS	86
APPENDIX C: CERTIFICATION AND COMPLIANCE	87
APPENDIX D: REPLACEABLE PARTS	88
INDEX	89

List of Tables

Table 2 - 1: TUNICS-Plus Specifications	15
Table 3 - 1: Models and ordering codes	19
Table 3 - 2: Options and ordering codes	20
Table 3 - 3: Accessories	20
Table 6 - 1: Crossover, 3-wire RS-232 serial connection to TUNICS- Plus.....	60
Table 6 - 2: Example of dialog between the computer and TUNICS-Plus (RS-232 C).....	60
Table 6 - 3: Examples of message formats with numerics (RS-232 C)	62
Table 6 - 4: Command parameters (RS-232 C).....	62
Table 6 - 5: Operating mode commands (RS-232 C).....	63
Table 6 - 6: System parameters (RS-232 C)	64
Table 6 - 7: Other system commands (RS-232 C)	65
Table 6 - 8: Supported GPIB functions.....	66
Table 6 - 9: Bit layout in status word (GPIB)	67
Table 6 - 10: Examples of single-indicator SRE values.....	68
Table 6 - 11: Sample command syntax (GPIB)	70
Table 6 - 12: IEEE Standard commands	72
Table 6 - 13: Parameter commands (GPIB)	72
Table 6 - 14: Operating mode commands (GPIB)	73
Table 6 - 15: System parameter queries (GPIB)	74
Table 7 - 1: Test equipment	77

List of Figures

Figure 2 - 1: TUNICS optical design	12
Figure 3- 1: Connect the power cord to the rear panel of TUNICS-Plus	21
Figure 4- 1: TUNICS-Plus front panel	27
Figure 4- 2: TUNICS-Plus rear panel.....	31
Figure 5 - 1: Instrument Setup for Wavelength Monitoring	46
Figure 5 - 2: Wavelength Monitoring Scheme.....	48
Figure 5 - 3: Auxiliary Input/Output Connectors on Rear Panel of TUNICS-Plus	50
Figure 5 - 4: Low frequency modulation setup	51
Figure 5 - 5: Resulting LF-modulated diode current.....	52
Figure 5 - 6: High frequency modulation setup	53
Figure 5 - 7: Resulting HF-modulated diode current	54
Figure 5 - 8: Remote Interlock SUBD-9 connector	56
Figure 7 - 1: Replacing the line fuse	82

ABOUT THIS MANUAL

Purpose

3

Audience

3

Prerequisites

3

Introduction

4

Chapters/Appendices

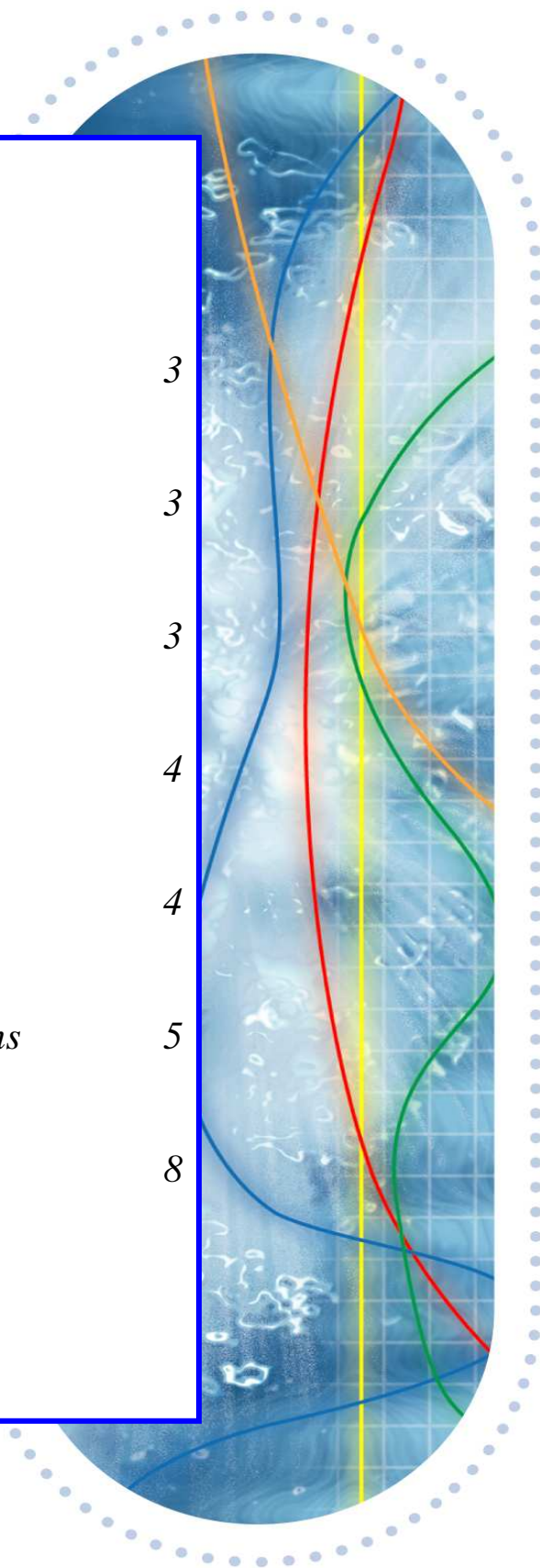
4

General Safety Considerations

5

Contacting Yenista

8



About this MANUAL

Purpose

The purpose of this instruction manual is to explain how to set up the TUNICS-**Plus** and to utilize the various features and the available control modes:

- TUNICS-**Plus** manual mode
- RS-232C and IEEE-488.1 remote modes

Information in this manual applies to TUNICS-**Plus** series including TUNICS-**Plus** O, **Plus** E, **Plus** S, and **Plus** CL models, as well as the wide-band (/WB) option.

Audience

This manual is for anyone wishing to use a TUNICS-**Plus** Tunable External Cavity Laser to test fiber-optic systems and components as well as all other performance-intensive laser source applications.

This manual's architecture and its contents are designed to appeal to professionals with the necessary technical background and prerequisites needed to use this type of laser equipment.

Prerequisites

For this instruction manual it is assumed that you:

- Are familiar with fiber optic technology.
- Are familiar with the RS-232 C and/or IEEE-488.1 interfaces used to control the TUNICS-**Plus** in the remote modes.

About this MANUAL

Introduction

This manual is divided into:

- Contents
- Chapters
- Appendices

To aid in searching for specific information, you will find clearly marked divider pages that indicate the various sections of the manual.

Chapters/Appendices



This document consists of six chapters designed to guide you through a correct installation/operation of your laser equipment. It also introduces theoretical/design considerations and provides servicing/maintenance information, as well. In addition, you will find four appendices that contain useful information on troubleshooting the system, the valid range of parameters setting, the product certification and declarations of compliance, and a list of replaceable parts.

Chapter/Appendix	Main Topics Covered
General Information	<ul style="list-style-type: none">• Product Overview• Optical Design• Benefits• Specifications
Getting Started	<ul style="list-style-type: none">• Unpacking• Standard Accessories• Options• Installation & First-Time Operation
Front & Rear Panels	<ul style="list-style-type: none">• Front Panel• Rear Panel
Operation	<ul style="list-style-type: none">• Introduction• Most Commonly Used Functions• Advanced Functions• Internal Wavelength Referencing• Active Cavity Control• Active Wavelength Monitoring• Auxiliary Inputs and Outputs
Remote Control	<ul style="list-style-type: none">• Introduction• RS-232 C Remote Control• IEEE-488.1 Remote Control
Service Information	<ul style="list-style-type: none">• Performance Verification• Maintenance
Appendices	<ul style="list-style-type: none">• Appendix A: Error and Warning Messages• Appendix B: Valid Range for Parameter Setting• Appendix C: Certification and Compliance• Appendix D: Replaceable Parts• Index

General Safety Considerations

Terms and Symbols in this Manual

The following terms and symbols may appear throughout this manual:

	WARNING	Warning statements identify conditions or practices that could result in injury or loss of life.
	CAUTION	Caution statements identify conditions or practices that could result in damage to this product or other property.







Terms on the Product

The following terms may appear on the product:

DANGER	Indicates an injury hazard immediately accessible as you read the marking.
WARNING	Indicates an injury hazard not immediately accessible as you read the marking.
CAUTION	Indicates a hazard to property, including the product itself.

Symbols and Labels on the Product

The following symbols and labels may appear on the product:

	CAUTION: refer to manual
	Risk of electric shock
	Protective Ground (Earth) Terminal
	AVOID EXPOSURE Label (located on the front panel.)
	LASER SYMBOL Label (located on the fiber-optic connector cover on the front panel.)
	Compound Label for Safety Label and Verification. Contains information such as Laser Class IIIb DANGER-Label, Electrical Warning, Fuse description, and Standard Compliance.

Laser Safety

The following tables provide information on general laser specifications, laser class, maximum wavelength ranges, and output power of the various TUNICS-**Plus** models and options.

General Laser Specifications

Laser type	External cavity laser diode
Laser output	Single mode optical fiber
Output divergence (full angle @ $1/e^2$)	0.18 rad

Laser Class

Standard	Laser Class
IEC 60825-1	Class 1M
21 CFR Subchapter J	Class IIb




Note: TUNICS-**Plus** is a Class IIb product (21 CFR Subchapter J) equipped with a Remote Interlock connector (See the *Auxiliary Inputs and Outputs* section). The laser source is classified as class 1M according to IEC 60825-1 (2001).

CLASS 1M LASER PRODUCT
(IEC 60825-1 / 2001)

Absolute Maximum Ratings

TUNICS Model	Plus O	Plus E	Plus S
Wavelength range (nm)	1240-1380	1330-1470	1410-1550
Max power rating	20 mW	20 mW	20 mW





TUNICS Model	Plus S WB	Plus CL	Plus CL WB
Wavelength range (nm)	1370-1560	1500-1640	1470-1660
Max power rating	20 mW	20 mW	20 mW

	WARNING	The use of controls or adjustments not specified herein as well as the performance of unauthorized procedures may result in hazardous radiation exposure.
	WARNING	The use of optical instruments with this product will increase eye hazard. Do not under any circumstances look directly into the fiber end of an optical cable attached to the optical output while the device is in use for this may cause permanent eye damage and possible loss of eyesight. Note that the laser radiation is not visible to the human eye, therefore, protective cap must always be replaced on the laser output connector after use to avoid involuntary exposure to laser radiation.
	WARNING	Disable the laser output before connecting or disconnecting a fiber optic cable on the instrument.


Equipment Safety Warning

If this equipment is mishandled, abused, or used in a manner not as specified by Yenista in this instruction manual, the protection provided by the equipment may be impaired.

Electrical Safety

	WARNING	<p>Power Supply Safety Requirements</p> <p>TUNICS-Plus is intended to operate from a power source that does not apply more than 265 volts RMS between the supply conductors or between either of the supply conductors and the ground.</p> <p>TUNICS-Plus has a chassis connected to ground via the power supply cable. A protective ground connection by way of the grounding conductor in the power cord is essential for safe operation.</p> <p>To avoid the possibility of injury, insert the power cable only into a socket outlet with a protective earth contact. Before switching on the instrument, make sure the electrical installation fulfills the local safety requirements.</p>
	WARNING	<p>Do not attempt to perform servicing or maintenance.</p> <p>To avoid personal injury, do not operate this equipment without the protective cover of the chassis. Do not make any service or maintenance of any kind to the system. Refer servicing to authorized Yenista personnel only.</p>
	WARNING	<p>Use the proper fuse.</p> <p>To avoid fire hazard, use only the correct type fuse, voltage, and current ratings as specified in the <i>Maintenance</i> section in the <i>Service Information</i> chapter.</p>
	WARNING	<p>Operate under proper environmental conditions.</p> <p>TUNICS is not designed for outdoor use. To avoid the possibility of injury, do not expose the instrument to rain or excessive moisture. Do not operate the instrument in the presence of flammable gases or fumes.</p>

Ventilation of the Unit

	CAUTION	<p>Make sure there is sufficient clearance below and at the back of the unit to ensure proper ventilation.</p> <p>The apertures (under the unit) and the fan (behind the unit) are used for the inside ventilation.</p> <p>When used as a stand-alone benchtop unit, this product is designed so that the feet should leave enough room under the unit to enable proper the ventilation.</p> <p>If you plan to rack-mount this unit, make sure to allow at least a 15-mm (0.6 inches) gap between two units.</p>
---	----------------	---

Contacting Yenista

Yenista Optics

4 rue louis de Broglie
BP 80429
22304 LANNION
France

Tel : +33 2 96 48 37 15

Fax: +33 2 96 48 73 04

E-mail: sales@yenista.com

www.yenista.com

GENERAL INFORMATION

Product Overview

11

Optical Design

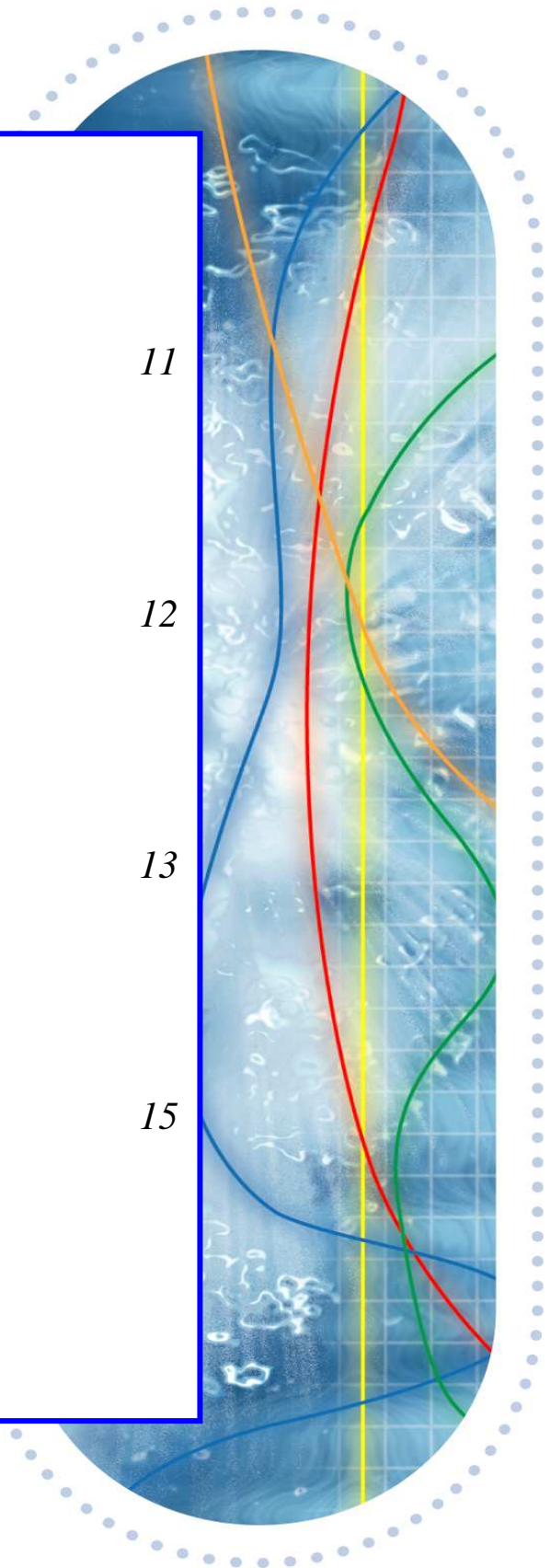
12

Benefits

13

Specifications

15



General Information

Product Overview

TUNICS-**Plus** is a high-performance, self-aligned cavity laser that implements Yenista state-of-the-art proprietary configuration. It is your ideal partner for all your design applications of your optical communications systems.

TUNICS-**Plus** is a general-purpose instrument designed for fiber-optic system testing and precise WDM component testing and for research and development in the field of optical communications. It offers maximum flexibility to cover all present and future requirements for tunable laser diode sources in the lab. TUNICS-**Plus** instruments are available in different versions providing a variety of wavelength and optical output power ranges. TUNICS-**Plus** is *the* fully integrated bench-top instrument fundamental to your lab applications. It comes with an optical head already integrated inside the electronics controller.

TUNICS-**Plus** comes as a complete line of tunable laser diode sources with wavelengths spectrum covering the S, C, and L bands, from as low as 1260 nm for your S-band applications up to 1640 nm for the L-band extended DWDM systems.

The sweeper mode enables TUNICS-**Plus** instruments to feature mode-hop-free, truly continuous wavelength sweep over the operational wavelength range.

The wide-band option (/WB, available upon order for S and CL models; see *Options and Ordering Information* on page 19) provides a 150 nm wide operational wavelength range¹.

TUNICS-**Plus** implements its own internal wavelength referencing system to provide ± 40 pm wavelength accuracy with a remarkable ± 5 pm repeatability. In addition, TUNICS-**Plus** tunable external cavity laser supports active wavelength monitoring via an external wavemeter. This setup provides outstanding accuracy to your tunable laser source within the accuracy of the wavemeter (virtually ± 2 pm wavelength accuracy, depending on the performance of the wavemeter). This makes the TUNICS-**Plus** laser source the ideal tool for precise WDM component characterization.

¹ @ 0 dBm output power

Optical Design

The optical layout is a modified Littman-Metcalf configuration. Double-pass reflection on the grating provides maximum dispersion, while the very short cavity, which maximizes mode spacing, yields a spectrally pure and truly single-mode operation.

The high-performance optical isolator and angle-polished output fiber connector also contribute to spectral purity by protecting the laser cavity from any interference that may be generated by your equipment.

In the TUNICS configuration, however, the dihedral rear reflector in the self-aligned cavity acts as a 2D *corner cube*, so that the resonator remains perfectly in tune, overriding any small misalignments and thus guaranteeing long-term stability.

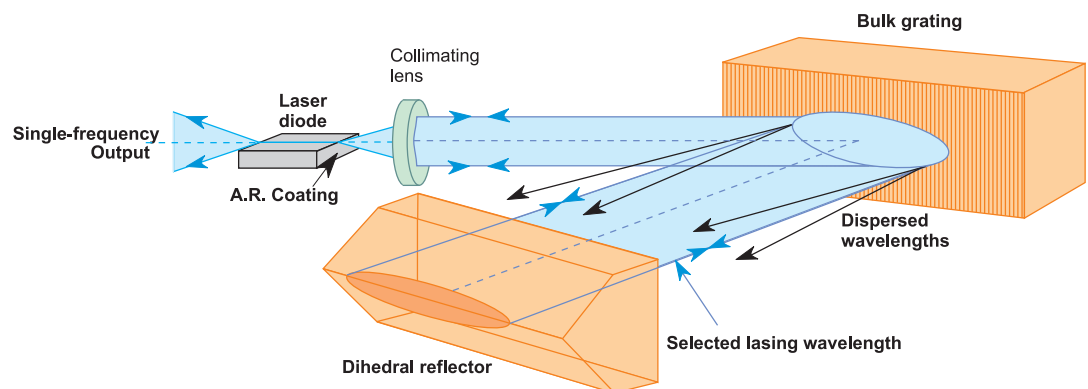


Figure 2 - 1: TUNICS optical design

Benefits

The proprietary optical design of the TUNICS-**Plus** Tunable External Cavity Laser provides a unique range of user benefits.

Benefit	Description
Wide, fast, and 100%-continuous tunability	Extremely smooth scans over 100 nm obtained in 1 second, with a high resolution of 1 picometer. A digitally-controlled analog fine-tuning extending the resolution beyond the one-picometer step to the sub-MHz level.
Active cavity control	A new proprietary active cavity control provides the means for full optical communications systems and components characterizations. It ensures true mode-hop-free operations over the entire tuning range along with remarkable wavelength accuracy.
Internal wavelength referencing	To ensure highly accurate (± 40 pm) and repeatable (± 5 pm, typ.) measurements, TUNICS- Plus features its own internal referencing system.
Active wavelength monitoring	Use the active wavelength monitoring feature, to turn your TUNICS- Plus into a high-precision tunable laser for the characterization DWDM components.
Full-featured wavelength sweeper mode	True swept-wavelength function with adjustable sweep speed setting over the user-defined wavelength range (within the system's valid wavelength range).
Mode-hop free	A wide range guaranteed to be free of any mode hopping to ensure a smooth and accurate wavelength sweep for reliable testing of narrow-band components.
High output power	More than 0 dBm guaranteed over the entire spectral range in all models. +8 dBm delivered by TUNICS- Plus S and CL over designated wavelength range (See <i>Specifications</i>).
Wide-band option (/WB)	TUNICS- Plus models featuring the wide band option can operate fully on a 150 nm extra-wide wavelength range (See <i>Specifications</i>).
Outstanding long-term stability	Self-aligned optical layout, single-moving-part design and all-invar construction to ensure a high level of long-term stability.
Intuitive and user-friendly controls	Optimized keyboard and display for intuitive and easy-to-implement laboratory operations. Dual parameter input (keyboard or multi-speed rotary control).

General Information

Benefit	Description
Multiple modulation capability	Easy analog and digital modulation of the optical power from DC to 1 GHz. Low frequency modulation from 30 kHz to 8 MHz. High frequency modulation from 30 kHz to 1 GHz. Mode-locked operation around 5 GHz. The 150 kHz linewidth can also be degraded to more than 100 MHz when high coherence is required.
Easy system integration	A fully integrated system featuring computer interfaces, onboard software, and analog inputs/outputs to allow all operations to be performed in remote mode and to make system integration easy.

Specifications

Table 2 - 1 lists the technical specifications of the various TUNICS-**Plus** Tunable External Cavity Laser.

Table 2 - 1: TUNICS-Plus Specifications

Tuning characteristics (1/2)	Plus O	Plus E	Plus S
Wavelength range (nm) (mode-hop-free)			
P = 0 dBm	1260-1340	1340-1430	1430-1530
P = 3 dBm	1280-1320	1380-1410	
P = 6 dBm			1440-1500
P = 8 dBm			1450-1490
Tuning characteristics (2/2)	Plus S WB	Plus CL	Plus CL WB
Wavelength range (nm) (mode-hop-free)			
P = 0 dBm	1390-1540	1525-1625	1490-1640
P = 3 dBm			
P = 6 dBm	1420-1520	1540-1620	1520-1630
P = 8 dBm	1440-1510	1560-1600	1540-1610
Step-mode characteristics			
Mode hop free range	Whole wavelength range for each specified power		
Absolute wavelength accuracy ²	±0.04 nm		
Tuning repeatability (typ.)	±0.005 nm		
Wavelength setting resolution	0.001 nm		
Optical frequency fine tuning	±2 GHz		
Tuning speed (typ.)	1 s (100 nm)		
Power stability ³	±0.01 dB		
Swept-mode characteristics (Optional /SW)			
Mode hop free range	Whole wavelength range for each specified power		
Scan speed	Adjustable from 1 to 100 nm/s		
Power flatness during scan (typ.)	±0.25 dB		
Power repeatability from scan to scan (typ.) ⁴	±0.05 dB		
Laser output characteristics			
Linewidth (coherence control OFF)	400 kHz typ.		
Linewidth (coherence control ON)	>100 MHz		
Side Mode Suppression Ratio ⁵	>45 dB (except Plus O & Plus E: > 40 dB)		
Signal to source spontaneous emission ratio ⁶	>55 dB (except Plus O & Plus E: > 45 dB)		
Relative Intensity Noise (RIN) ⁷	-145 dB/Hz (typ.)		

Specifications are given after 2 hours warm up.

² After self calibration.

³ Over one hour at constant temperature.

⁴ Over 100 scans at constant temperature.

⁵ Measured with 0 dBm output power.

⁶ Spontaneous emission measured on a 0.1 nm bandwidth at ±1 nm from the signal.

⁷ Measured at an electrical frequency of 100 MHz.

Table 2 - 1: TUNICS-Plus Specifications (cont'd)

Interface	
Optical interface	FC-APC connector
Output fiber	SMF-28™
Output isolation	35 dB
Return loss	60 dB
Remote control	IEEE-488-1 and RS-232 C
Low frequency modulation	
• APC mode	10 kHz to 8 MHz
• Direct current mode	DC to 8 MHz
High frequency modulation	30 kHz to 200 MHz
Mode-lock frequency	5 GHz
Environment	
Operating temperature range	+15 to +30 °C (+60 to +85 °F)
Power supply (auto-switch)	100 to 240 V, 50 to 60 Hz
Dimensions (W x H x D)	448 x 133 x 370 mm ³
Weight	12.5 kg


GETTING STARTED

Unpacking 19

Standard Accessories 19

Options 19

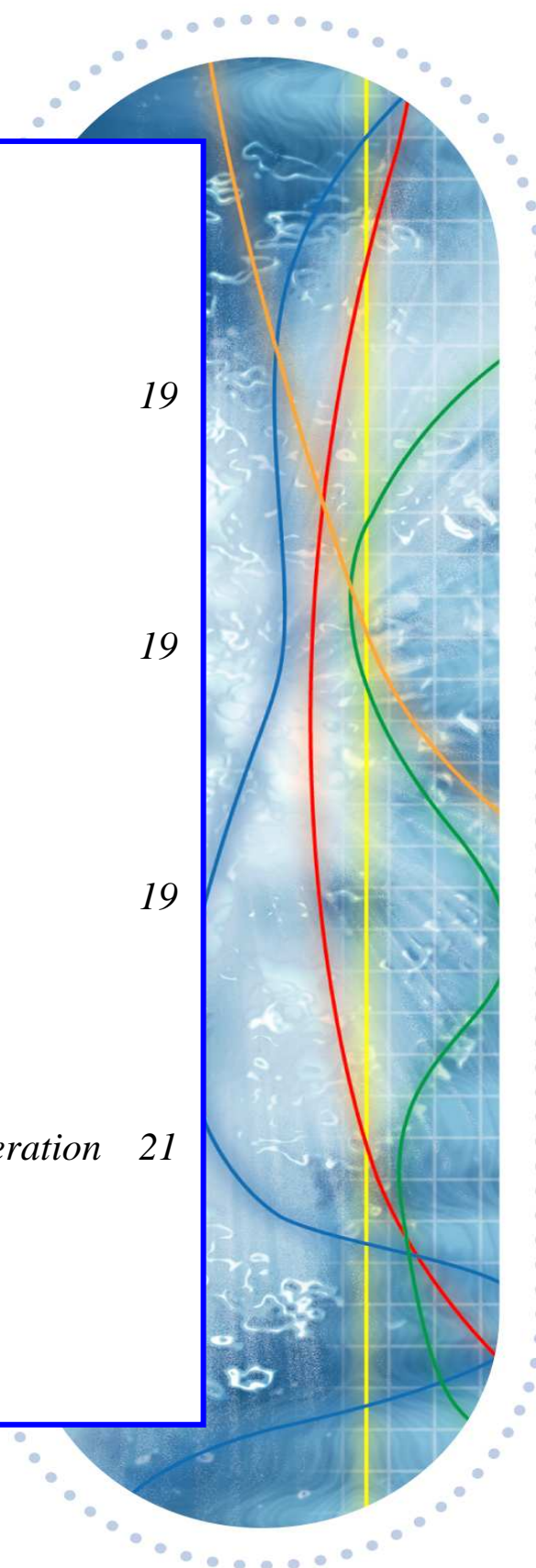Installation & First-Time Operation 21



Getting Started

Unpacking

Handle the device with care when unpacking. To pull out the unit from its packaging, hold it by its two lateral retractable handles to help you keep the unit horizontal. After unpacking, set the device on a flat stable surface that is free of excessive vibration.



CAUTION

Keep the original shipping container for use in case the instrument needs to be returned to Yenista for repair or servicing.

Standard Accessories

The TUNICS-Plus Tunable External Cavity Laser comes with the following standard accessories:

- Certificate of traceable calibration
- Acceptance test report (ATR)
- FC/PC to FC/APC adapter cable
- Safety-lock key
- Instruction manual
- Power cord
- Rack-mount handles (2x)

Options and Ordering Information

TUNICS-Plus comes in various models and options available upon order by specifying the corresponding ordering code. Accessories may be purchased separately.

Models

Model names and main characteristics are listed below:

Table 3 - 1: Models and ordering codes

TUNICS-Plus Model Ordering Code	Model characteristics
O	O-band over 1260-1340 nm; +3 dBm output power.
E	E-band over 1340-1430 nm; +3 dBm output power.
S	S-band over 1430-1530 nm; +8 dBm output power.
S WB	S-band over 150 nm wide-band wavelength range of 1390-1540 nm; +8 dBm output power.
CL	C+L-band over 1525-1625 nm; +8 dBm output power.
CL WB	C+L-band over 150 nm wide-band wavelength range of 1490-1640 nm; +8 dBm output power.

Options

Available options include:

Table 3 - 2: Options and ordering codes

Option Ordering Code	Option
WB	Wide-band ⁸ 150 nm wavelength range @ 0 dBm.
M	Polarization maintaining output fiber (orientation TE in slow axis, in line with connector key)

Note: For the standard benchtop model, you do not need to specify any option.

Accessories

Available accessories include:

Table 3 - 3: Accessories

Labview driver for TUNICS- Plus
Fiber-optic jumper FC-APC/FC-APC
Fiber-optic jumper FC-APC/FC-PC (delivered as standard with each instrument)
Polarization maintaining fiber-optic jumper FC-APC/FC-APC
Polarization maintaining fiber-optic jumper FC-APC/FC-PC
Field carrying case (3U)

Ordering code

Ordering information must be provided in the following manner, as shown in the table below (as an example only):

Instrument name	Family	Model	Option 1	Option 2
TUNICS	Plus	/ S	/ WB	/ M

⁸ WB option available for TUNICS-**Plus** S and CL only.


Installation & First-Time Operation

To remove the unit from place to place, use the handles on both sides which help to keep the instrument horizontal.

Installation

To install TUNICS-**Plus** tunable laser source proceed as follows:

1. Set the laser source on a flat stable surface that is free of excessive vibration.
2. Allow the flow of air from the cooling fan to circulate freely around the device and remove any equipment or paper that could block the air-flow.
3. Use the proper power cord for your location.
4. Make sure the key-lock switch on the front panel is in the Standby mode before applying power.
5. Connect one end of the power supply cable to the rear panel of the unit and plug the other end to the proper voltage mains supply point. TUNICS-**Plus** is equipped with a self-regulating power supply that adapts to both AC 110V and 230V voltages and 50 Hz or 60 Hz frequencies.
6. To turn ON the unit, set the main switch located on the rear panel to the "I" position.

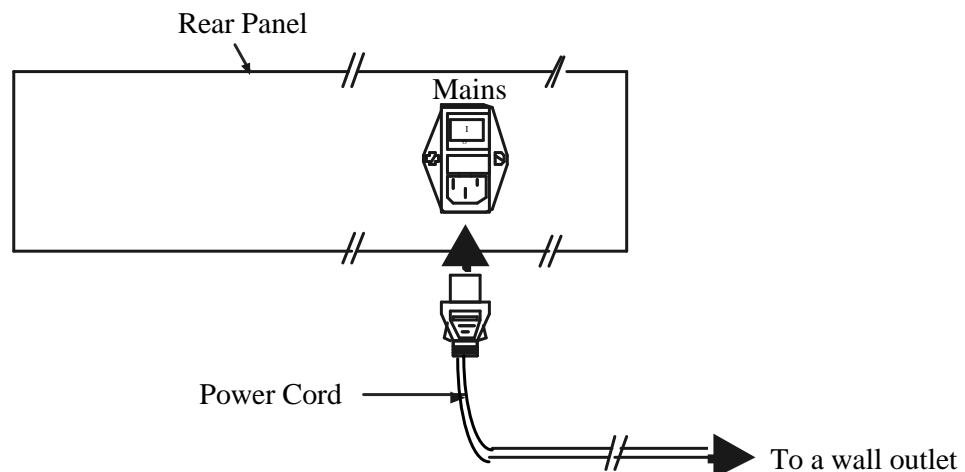

*Figure 3- 1: Connect the power cord to the rear panel of TUNICS-**Plus***

Initializing the Laser

For safe operation, use the following procedure to initialize the TUNICS-**Plus** laser:

1. Make sure the key-lock switch on the front panel is in the Standby position.
2. Set the rear panel MAINS switch to ON (I) (see Figure 3- 1).
3. Turn the front panel key to the right and set it to ON. During the initialization routine, the message “**Initializing...**” is displayed.
4. The message “**Referencing...**” is then displayed to indicate that the system performs its self-calibration (internal wavelength referencing). This message clears once the laser is initialized and ready for operation.
5. Once the initialization routine is complete, the laser output is ready to enable. The message “**Disabled**” is displayed. The cavity laser is tuned on its central emission wavelength (depending on the model).

TUNICS-**Plus** performs self-checks during initialization sequence. The APC LED is lit. Power setting is set to zero when the laser is powered up.

Note: User can decide to skip the self-calibration (skip Step 3 - Referencing) procedure by simply pressing the “0” key while the message “Initializing...” is displayed. Note that for optimum system performance, we recommend you to let TUNICS-**Plus** perform its self-calibration routine when turning on the system.

	CAUTION	Avoid switching the unit off before the system is fully initialized.
--	----------------	--

Enabling the Laser Output

For laser safety reasons, the TUNICS-**Plus** laser output is not immediately enabled after system turn-on and initialization. After the initialization routine, TUNICS-**Plus** is in Disabled mode. TUNICS-**Plus** features an Enable key to enable/disable laser output.

- To enable the laser output, press the **Enable** button. The Enable key is lit.
- To disable the laser output, press the **Enable** button again and check that the Enable key LED is off.

Optical Connection to the Unit

External optical equipment can connect to TUNICS-**Plus** via the FC/APC fiber-optic connector located on the front panel in the Input area.

	CAUTION	Keeping fiber-optic connectors clean at all times is essential to achieve optimum system performance. Refer to the cleaning instructions in the <i>Maintenance</i> section. The standard output connector is a FC/APC type connector. Never connect another type of connector to the optical output.
--	----------------	--

Shutting Down the TUNICS-Plus Laser

To shut down the unit, follow these steps:

1. Press the **Enable** button to disable the laser output and make sure the Enable key LED is off.
2. Turn the front panel key to the Standby position (counterclockwise).
3. The message “Parking...” is displayed and the optical head is automatically moved to its parking position.
4. Wait a few seconds until complete system shutdown.

	CAUTION	During system shutdown, DO NOT TURN OFF the system with the rear panel O/I switch.
	CAUTION	Before moving the unit, or when the unit is not in use for a long period of time, the rear panel Mains O/I switch must be set prior to the “O” position. Be sure the system is completely turned off before transporting it
	CAUTION	If the TUNICS-Plus unit has been turned off, you must wait at least thirty seconds before initializing the laser again.

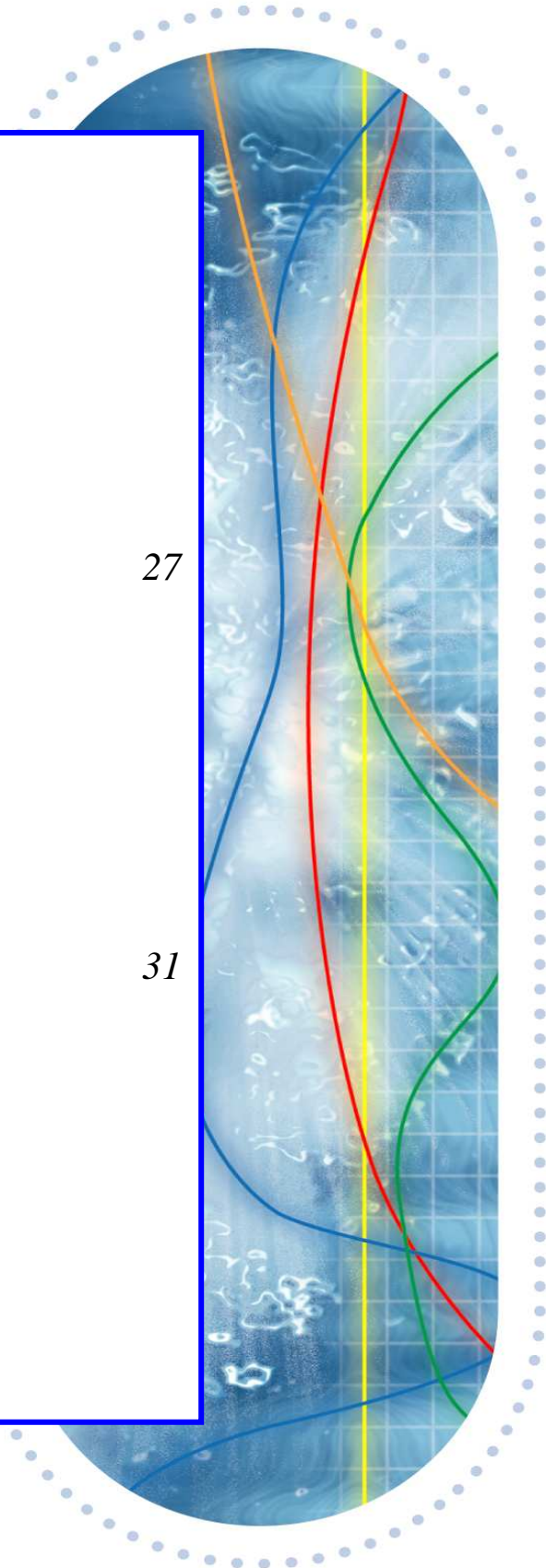
FRONT & REAR PANELS

Front Panel

27

Rear Panel

31



Front & Rear Panels

This section provides an overview of the operating controls of all TUNICS-Plus functions.

Front Panel

Figure 4- 1 shows the front panel of the TUNICS-Plus Tunable External Cavity Laser.



Figure 4- 1: TUNICS-Plus front panel

Description of the Front Panel Functions Keys

The TUNICS-Plus front panel is divided into two areas, as illustrated in Figure 4- 1:

DATA area	To display the wavelength, optical power, diode current, and messages.
MODE area	To select the operating mode. To select the system parameter to be modified.
SETUP area	User interface to enter the parameter value.

The various function keys are described below:



DATA Area

The DATA area displays the two main operating parameters, the emission wavelength (λ) (or optical frequency f), the output power (P), and the laser diode current (I); it is also used to display system-status messages.

Function Key	Purpose
nm/GHz	Displays either the wavelength in nanometer (nm) or the optical frequency in GigaHertz (GHz).
mW/dBm	Displays the power either in dBm or mW.
Enable	Enables the optical output. (Note that in Disable mode, the message “Disabled...” is displayed.)

MODE Area




The **MODE** area includes 8 function keys that are lit when the corresponding function is active. These function keys allow the user to access and change parameters that need be configured (see the *General Parameter Setting Principle* section that follows).

λ, f	Used to edit the current value of either the emission wavelength λ or the optical frequency f , and to allow entry of a new value.
P	Used to edit the current value of the power output P , and to allow entry of a new value. Switches to the Constant-Power mode.
I	Used to edit the current value of the diode current I , and to allow entry of a new value. Switches to the Constant-Current mode.
STEP	Used to change the wavelength increment (STEP) value. It is also possible to increase or decrease the wavelength STEP-by-STEP, at any time, by pressing the  or  arrows on the keypad.
APC	(Automatic Power Control) switches between constant current (APC light off) and constant power (APC light on) modes.
FSC	FSC (Fine SCanning mode) enables a complete and virtually continuous sweep of the wavelength on a ± 2 GHz range. Once this mode is selected, the wavelength value can be adjusted by rotating the Modify knob located in the SETUP area. The wavelength sweep is displayed in the DATA area either in nm or in GHz accordingly. To exit the fine scanning mode, press the FSC key again
Remote	Press this button to switch between remote and manual modes when the system is connected to a remote terminal or a computer via the RS-232 C or IEEE-488.1 interface. The Remote button LED is lit when data is received by TUNICS-Plus on one of its RS-232 C or IEEE-488.1 ports. When the Remote key is on, the front panel keys become inactive. To regain manual control of the TUNICS-Plus at the front panel, press the Remote button again (the Remote LED is off).
2nd	The 2nd function key allows access to advanced functions (see the <i>Advanced Functions</i> section in the <i>Operation</i> chapter).
2nd + λ, f	To set the motor-speed for swept-mode wavelength scanning.
2nd + P	To perform power calibration.
2nd + I	To perform system auto-calibration.
2nd + APC	To disable/enable active cavity control.
2nd + dBm/mW	To turn on active wavelength monitoring.
2nd + STEP	To activate wavelength step-mode or sweep-mode scanning
2nd + 1	To disable/enable backlash suppression of the micrometer screw.
2nd + 2	To display the software release.
2nd + Remote	To change the GPIB address. Also to activate/deactivate the RS-232 remote mode. Note: to reach this option, you first need to confirm the GPIB address (via the Enter key).

2nd + FSC	<p>Enables the coherence control function. The Coherence Control function adds noise to the optical frequency output of the laser producing an apparent broadening of the line-width. Coherence control is particularly useful to eliminate parasitic interferometers in the user setup. Note: This function can also be accessed by using the FSC input on rear panel.</p> <p>(See <i>Auxiliary Inputs and Outputs</i> section in <i>Operation</i> chapter.)</p>
------------------------------------	--

SETUP Area

The SETUP area is used to enter or change the values of the system operating parameters. It includes:

- Numeric keypad
-  arrow key to correct a typing error
- **Enter** key to validate the choice
-  key to enter a decimal point
-  to enter a negative value
- Rotating knob to adjust parameter values

In order to change a parameter, first select the appropriate parameter key in the MODE area. The display shows the current setting for the selected parameter, which can be readily modified by entering the new value directly at the numeric keypad or by adjusting the parameter to the desired value with the rotating Modify knob.



General Parameter Setting Principle

The various parameters may be set or changed by using one of two methods:

1. At the numeric keypad
2. By using the rotating Modify knob



Using the Numeric Keypad

To set a parameter using the numeric keypad, follow these steps:




Step	Action
1. <u>Result:</u>	Press the appropriate function key in the MODE area. The corresponding LED indicator is lit. The current parameter value is displayed. The = sign flashes to indicate that the value can be changed.
2.	<p>Enter the required value via the keypad.</p> <p>If you make a mistake while entering the value, use the arrow keys  /  to correct your entry.</p> <p>If the value entered is outside the valid range for the parameter selected, the following error message is displayed: Value error, press Ent Press Enter to continue. Note: The value is not changed and any scanning or calibration operation in progress is suspended.</p>
3.	To validate your change, press the same function key in the MODE area again.
4.	Apply the new parameter value by pressing the Enter key.

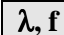Using the Rotating Modify Knob


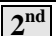

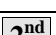

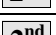
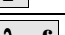


By rotating the Modify knob, you can update a parameter value in real time.
To set a parameter using the Modify knob, follow these steps:

Step	Action
1. <u>Result:</u>	Press the appropriate function key in the MODE area. The corresponding LED indicator is lit. The current parameter value is displayed. The = sign flashes to indicate that the value can be changed.
2.	Rotate the Modify knob until the required value is displayed.
3. <u>Result:</u>	To change the basic increment/decrement applied by the Modify knob, press the  and  arrows. A similar message is displayed: Scroll 0.01 Indicates which digit can be changed using the Modify knob.
4. <u>Result:</u>	To validate your change, press the same function key in the MODE area again. The new value is applied automatically. Important: The backlash suppression feature applied to the tuning micrometer screw is not activated when you change a value via the Modify knob, as the parameter is changed in real time. This is indicated by a dot at the end of the wavelength displayed similar to the following. Example: $\lambda = 1555.000 \text{ nm} \bullet$

The numeric keypad or control knob can edit the following parameters:

	Emission wavelength
	Optical power
	The diode current

User may also edit other parameters in the following manner:

	Wavelength step
 + 	To activate wavelength step-mode or sweep-mode scanning
 + 	Power calibration
 + 	To set the speed of the wavelength swept-mode scanning
 + 	To change the GPIB address Note: also used to activate/deactivate the RS-232 remote control mode.

Rear Panel

Figure 4- 2 shows the rear panel of the TUNICS-Plus.



Figure 4- 2: TUNICS-Plus rear panel

Description of the Rear Panel Connectors

The MAINS module includes the mains power switch, power cord connection, and fuse drawer.

In addition to the MAINS module, the following connectors are available on the rear panel:

Monitor (Scaled) Outputs

- Sync (motor Sync)
- λ (wavelength)
- P (power)
- I (current)

TTL Outputs

- Mod Sync. modulation sync (TTL output)
- Mod Out modulation output (TTL output)

Inputs

- HF Mod. high-frequency modulation input (SMA connector)
- LF Mod. low-frequency modulation input (BNC connector)
- FSC (fine scanning and coherence control)
- Remote Interlock connector as required for Class IIb laser products


Remote Interfaces

- RS-232 C connector for the remote control mode (SUBD-9).
- IEEE-488.1 connector for the remote control mode.

For a detailed discussion on particular connectors, see the *Auxiliary Inputs and Outputs* section in the *Operation* chapter.

OPERATION

<i>Introduction</i>	35
<i>Most Commonly Used Functions</i>	36
<i>Advanced Functions</i>	38
<i>Internal Wavelength Referencing</i>	43
<i>Active Cavity Control</i>	45
<i>Active Wavelength Monitoring</i>	46
<i>Auxiliary Inputs and Outputs</i>	50

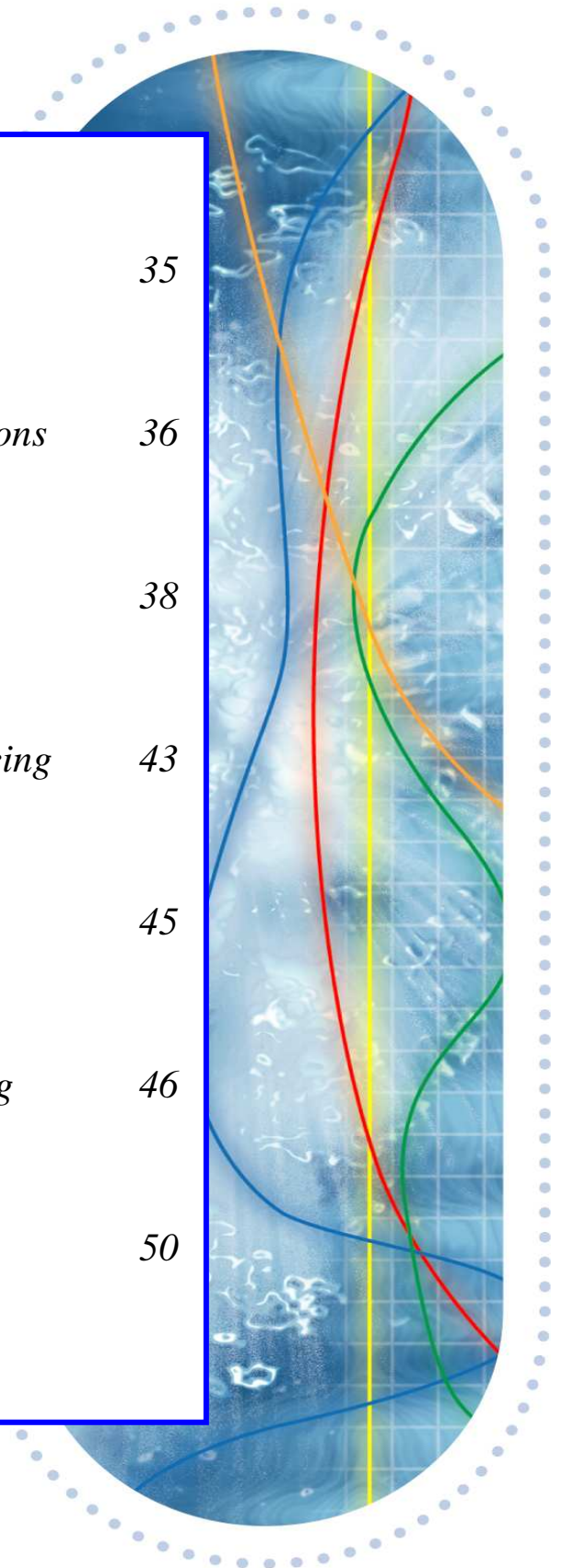

Operation

Introduction

This section describes how to use TUNICS-**Plus** Tunable External Cavity Laser in manual mode. TUNICS-**Plus** front panel allows direct access to most commonly used functions by pressing the corresponding function keys in the MODE area. In addition, a panel of advanced functions can be accessed via the 2nd function key.

The most commonly used functions are described as follows:

- To change the emission wavelength
- To change the wavelength step
- To change the beam's optical output power
- To change the diode current
- To switch between the constant-power and constant-current modes
- To use the fine scanning mode

The following advanced functions are also available:

- To implement step-mode wavelength scanning
- To implement swept-mode wavelength scanning (hop-free)
- To perform power calibration
- To turn on backlash suppression for the micrometer screw
- To change the GPIB address
- To turn on coherence control

TUNICS-**Plus** brings in many technological enhancements and new features. TUNICS-**Plus** enhanced features include:

- Internal wavelength referencing
- Active cavity control
- Active wavelength monitoring

Finally, auxiliary input and output are described at the end of this section.

Most Commonly Used Functions

To Change the Emission Wavelength

To change the emission wavelength, proceed as follows:

Step	Action
1.	Press the λ, f function key in the MODE area.
2.	Enter the new value using the numeric keypad or modify knob. Make sure to enter values within the valid range for parameter setting depending on your particular TUNICS-Plus model, as specified in <i>Appendix B: Valid Range for Parameter Settings</i> .
3.	Press the Enter key in the SETUP area.
	You may also increase or decrease the wavelength by pressing the ← / → arrow keys in the SETUP area, which has the same effect as rotating the modify knob.

To Change the Wavelength Step

To change the wavelength step, proceed as follows:

Step	Action
1.	Press the STEP function key in the MODE area.
2.	Enter the new value using the numeric keypad or modify knob. Make sure to enter values within the valid range for parameter setting depending on your particular TUNICS-Plus model, as specified in <i>Appendix B: Valid Range for Parameter Settings</i> . Note: It is possible that the required output power cannot be reached even at the maximum allowable current. If this occurs, TUNICS-Plus automatically limits the current to its maximum value and the message “Lim” appears at the right end of the displayed power, as shown in this example: P = 5.4 mW Lim
3.	Press the Enter key in the SETUP area.

To Change the Optical Output Power

To change the optical output power, proceed as follows:

Step	Action
1.	Press the P function key in the MODE area.
2.	Enter the new value using the numeric keypad or modify knob. Possible output power values range within 0.2 to 10 mW for the standard-power TUNICS-Plus model. Depending on which output power model you are using, make sure to enter values within the valid range for parameter setting of your particular TUNICS-Plus model, as specified in <i>Appendix B: Valid Range for Parameter Settings</i> . Note: It is possible that the required output power cannot be reached even at the maximum allowable current. If this occurs, TUNICS-Plus automatically limits the current to its maximum value and the message “Lim” appears at the right end of the displayed power, as shown in this example: P = 5.4 mW Lim
3.	Press the Enter key in the SETUP area.

To Change the Diode Current

To change the diode current, proceed as follows:

Step	Action
1.	Press the I function key in the MODE area.
2.	Enter the new value using the numeric keypad or modify knob. For possible values of the diode current, please refer to the test report provided with the system.
3.	Press the Enter key in the SETUP area.

To Switch Between Constant-Power and Constant-Current Modes

Pressing the **APC** function key allows the user to toggle between the constant-current and the constant-power modes.

If the APC key LED is...	Then the active mode is in...
On	<p>Constant-power mode.</p> <p><u>Result:</u> The laser-diode current is controlled to ensure an optical power output equal to the P parameter. By pressing the P function key, the system switches automatically to the Constant-Power mode. In this mode, it may not be possible to obtain the required power, even with the maximum allowable current. In this case:</p> <ol style="list-style-type: none"> 1. The current is limited to its maximum value. 2. The optical power is lower than the user-set power value. 3. The Lim message is displayed to the right of the current and power values.
Off	<p>Constant-current mode.</p> <p><u>Result:</u> The diode current level is determined by the I parameter. By pressing the I function key, the system switches automatically to the Constant-Current mode. Note: Using the P or I keys switches the system automatically in the Constant-Power or Constant-Current mode, respectively.</p>

To Use the Fine Scanning Mode


To...	Press the...
Enable the Fine Scanning mode	Press the FSC function key.
<u>Result:</u>	<p>The wavelength can be adjusted at all times using the rotating Modify knob.</p> <p>The wavelength sweep is displayed in the DATA area in nm or GHz, as required.</p>
Exit the Fine Scanning mode	Press the FSC function key again.

Advanced Functions

To Implement Step-by-Step Wavelength Scanning

Wavelength scanning enables the TUNICS-**Plus** tunable laser to broadcast on a different wavelength at each time interval. It scans incrementally the entire wavelength scanning range specified by the user.

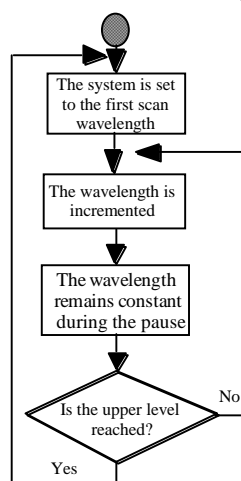The wavelength scan process is illustrated in the following flowchart:



Starting from the first wavelength (lower wavelength limit set by the user), the emission wavelength increases incrementally (the wavelength increment is stored in the Step parameter and set by the user) up to the last wavelength value (upper wavelength limit set by the user). During each scan interval, the TUNICS-**Plus** laser broadcast at a set wavelength for the duration specified by the Pause Timer parameter (Pause Timer set by the user).

Before activating the wavelength scanning, user must set the optical output power of the TUNICS-**Plus** laser (see the *To Change the Optical Output Power* section in this chapter).

To implement wavelength scanning, follow these steps:

Step	Action
1.	Press 2nd + STEP to activate the scan mode.
<u>Result:</u>	You are prompted to enter the lower wavelength scan limit while the current setting is displayed.
2.	Enter the lower wavelength scan limit and press Enter .
<u>Result:</u>	You are prompted to enter the upper wavelength scan limit while the current setting is displayed.
3.	Enter the upper wavelength scan limit and press Enter . (For valid scan limits, see Appendix B.)
<u>Result:</u>	You are prompted to enter the scan STEP while the current setting is displayed.
4.	Enter the scan STEP and press Enter . (To know the valid range for scan STEP values, see Appendix B.)
<u>Result:</u>	You are prompted to enter the Pause Timer between two scan steps while the current setting is displayed.

Step	Action
5.	Enter the Pause Timer and press Enter . Possible values for Pause Timer range from 0.1 to 25 seconds.
<u>Result:</u>	Scanning begins.
6.	To suspend the current wavelength scan, press the STEP key. You are then prompted to “stop scan 1/0?”. To stop the scan press the numeric keys 1 + Enter . Note: Scanning stops and the emission wavelength remains set to the value applied at the time the scan was suspended. To resume scanning, press numeric keys 0 + Enter .

In manual operating mode, the wavelength scanning process repeats continuously until user presses the **STEP** key to suspend or stop the scan.

⇒ Note that operating parameters cannot be changed while the wavelength scanning is in progress. Therefore, make sure the power output level you set can be achieved within the maximum current limits over the entire scanning range. If the desired power output cannot be reached due to current limitations, the **Lim** message is displayed to indicate the system no longer operates at constant power.

To Implement Continuous Wavelength Sweep (hop-free)

TUNICS-**Plus** allows you to perform truly continuous, mode-hop-free scans in swept-mode over the entire operational wavelength range.

Wavelength sweep commands on front-panel

To implement wavelength sweep, follow these steps:

Step	Action
1.	Press 2nd + STEP to activate the sweep mode.
<u>Result:</u>	You are prompted to enter the lower wavelength scan limit while the current setting is displayed.
2.	Enter the lower wavelength scan limit and press Enter .
<u>Result:</u>	You are prompted to enter the upper wavelength scan limit while the current setting is displayed.
3.	Enter the upper wavelength scan limit and press Enter . (For valid scan limits, see Appendix B).
<u>Result:</u>	You are prompted to enter the scan STEP while the current setting is displayed.
4.	When prompted, enter “0” as scan STEP parameter, and press Enter.
5.	When prompted, enter the PAUSE TIMER parameter, and press Enter. The PAUSE TIMER is the dwell time at the upper wavelength of the sweep range before the TUNICS- Plus tunable laser returns to the lower wavelength to restart a new sweep cycle. Possible values for Pause Timer range from 0.1 to 25 seconds.
<u>Result:</u>	The wavelength sweep (hop-free) is activated. Sweep begins.
	To suspend the current wavelength sweep, press the STEP key. You are then prompted to “stop scan 1/0?”. To stop the scan press the numeric keys 1 + Enter .

	<p>Note: Sweep stops and the emission wavelength remains set to the value applied at the time the scan was suspended.</p> <p>To resume scanning, press numeric keys 0 + Enter.</p>
--	---

In manual operating mode, the wavelength sweeping process repeats continuously until user presses the **STEP** key to suspend or stop the sweep. However, in remote operating mode (see *Remote Control* chapter), the system performs a single sweep path over the wavelength sweeping range and stops when the wavelength reaches the upper wavelength limit.

Sweep speed

You can set the sweeping process to scan from a lower to an upper wavelength value, at a user-specified speed. When the TUNICS-**Plus** has reached the upper wavelength limit, it returns to the initial wavelength value at maximum speed (i.e. 100 nm/s), and starts again the scanning process at the defined speed.

To set the sweep speed, follow these steps:


Step	Action
1.	Press 2nd + λ, f to display the current sweep speed setting.
2.	Enter the new sweep speed directly in “nm/s” unit and press Enter to confirm the change.

The sweeping speed is user-selectable within the speed range of 1 to 100 nm/s. You may choose any value of speed within this range. However, operational speeds that the motor can move are: [1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 17, 18, 20, 22, 25, 29, 33, 40, 50, 67, 100] (nm/s). Therefore, the speed you select is rounded-up to nearest value of operational speed.

For example, if you select 70 nm/s sweep speed at the TUNICS-**Plus** front panel, your setting is automatically corrected by the system to match the nearest operational value: 67 nm/s.

To Perform Power Calibration

The power calibration procedure calibrates the laser output power. The procedure requires the use of an optical power meter.

	<p>CAUTION</p> <p>This procedure clears the factory-set calibration data. Do not perform this procedure unless you actually intend to erase this data.</p>
---	---

To perform power calibration, follow these steps:

Step	Action
1.	Connect the laser output of the tunable laser to a power meter, and set the Enable mode (laser on).
2.	Press the 2nd + P function keys to activate the power calibration procedure.
<u>Result:</u>	The tunable laser moves to a short wavelength within the tuning range and displays the power level corresponding to the current calibration data.

3.		Measure the power level at the first wavelength on the power meter.
4.	<u>Result:</u>	Enter the reading at the prompt (P =). The tunable laser moves to a long wavelength within the tuning range and displays the power level corresponding to the current calibration data.
5.		Measure the power level at the second wavelength.
6.	<u>Result:</u>	Enter the reading at the prompt (P =). The power calibration is now complete and the instructions reference the new power calibration data.

To Turn On Backlash Suppression of the Micrometer Screw

The TUNICS-Plus firmware implements systematically an onboard procedure to eliminate the effect of backlash on the tuning micrometer screw. Therefore, you do not need to activate backlash suppression, which is already ON.

To Disable/Enable backlash suppression, follow these steps:

Step	Action
To disable the backlash suppression mode <u>Result:</u>	Press the 2nd + 1 function keys. A dot is displayed to the right of the wavelength. Example: $\lambda = 1555.000 \text{ nm}.$
To enable the backlash suppression mode	Press the 2nd + 1 function keys again.

To Change the GPIB Address

The default value for the GPIB address is 10.

To change the GPIB address, follow these steps:

Step	Action
1. <u>Result:</u>	Press the 2nd + Remote function keys. The current GPIB address is displayed as you are prompted to enter the new value at the numeric keypad.
2.	Enter the new GPIB address.
3. <u>Note:</u>	Press Enter to confirm the change. You are then prompted to activate/deactivate the RS-232 remote mode. Press Esc or refer below for further information.

To Turn On/Off the RS-232 Remote Control Mode

To activate/deactivate the RS-232 remote control mode:

Step	Action
1.	After confirming the GPIB address (refer to previous paragraph), you reach the status of the RS-232 remote control mode.
2.	You are then prompted to specify the mode status at the numeric keypad: <ul style="list-style-type: none"> Press “0” to turn the RS-232 mode “OFF”. Press “1” to turn the RS-232 mode “ON”.
3.	Press Enter to confirm the change.

To Turn On Coherence Control


To activate the Coherence Control feature:

Step	Action
1.	Press the 2nd + FSC function key. The FSC key LED then blinks. Coherence control provides an apparent broadening of the laser linewidth, which is useful in clearing out parasitic interferometers in the user setup.

Internal Wavelength Referencing

The TUNICS-**Plus** internal referencing system ensures optimum calibration of the laser source as early as system start-up.

Note: We also recommend user to perform an auto-calibration after the first half-hour of system warm-up or use.

	IMPORTANT	Change of Operating Environment When the system is subject to substantial changes of operating environmental conditions such as a sharp change in room temperature or mechanical vibrations, user must run a new system auto-calibration procedure, as specified below.
---	------------------	---

Self-calibration at start-up

TUNICS-**Plus** automatically performs a self-calibration procedure (also called Internal Referencing) during the initialization routine, which unfolds as follows:

Step	Action
1.	Turn the front-panel key to the right to begin system initialization.
2.	The system starts the initialization routine and the message “ Initializing... ” is displayed.
3.	The message “ Referencing... ” is then displayed to indicate that the system performs its self-calibration (internal referencing of the wavelength).
4.	After the initialization routine is complete, the message “ Disabled... ” is displayed which indicates that for safety reasons although the system is ready to use, the laser output is still off.

Note: User can decide to skip the self-calibration (skip Step 3 - Referencing) procedure by simply pressing the “0” key while the message “Initializing...” is displayed. Note that for optimum system performance, we recommend you to let TUNICS-**Plus** perform its self-calibration routine when turning on the system.



User-command for auto-calibration

In addition to the self-calibration performed at start-up only, user can instruct the system to perform a new auto-calibration at any time.

Follow these steps to perform system auto-calibration:

Step	Action
1.	Press the key the 2nd + I combination to start auto-calibration.
2.	TUNICS- Plus performs its internal referencing routine and the message “ Referencing... ” is displayed.
3.	After calibration is complete, TUNICS- Plus restores the initial wavelength value in effect at the time auto-calibration was performed.

Note: If auto-calibration has failed, the following message is displayed: **“Referencing error. Press Enter”**. Press Enter key and calibration is canceled. The system's current internal wavelength referencing remains unchanged.

	<p>IMPORTANT</p>	<p>Laser Safety & Auto-calibration</p> <p>While the system is performing an auto-calibration, the TUNICS produces a laser beam with an output power of 1 mW, even when the system is in <i>Disable</i> mode.</p> <p>To avoid eye damage, do not look into the laser output or into the end of an optical cable. Even though the laser radiation is not visible, the intense infrared light can cause eye damage.</p>
	<p>IMPORTANT</p>	<p>Wavelength Referencing & Remote Interlock</p> <p>For laser safety reasons TUNICS-Plus features a Remote Interlock connector at the rear-panel that can connect to an external remote interlock switch. When the interlock switch is open, the laser output is turned off. To enable the laser output, the interlock switch must be closed and user must press the Enable key.</p> <p>Remote interlock sets the system current and power to 0 mA and 0 mW, respectively. Therefore, while in remote interlock mode, the system cannot perform any self-calibration or internal referencing.</p> <p>When the system resumes from remote interlock mode into active state, we suggest user to run a new system auto-calibration procedure since the present system calibration data is still based on the previous system operating conditions.</p>

Active Cavity Control

TUNICS-**Plus** is equipped with onboard active cavity control.

When you turn on the TUNICS-**Plus**, after the initialization routine is complete, the system switches automatically to APC (Constant Power) mode with active cavity control (the APC LED is plain on).

Active cavity control operating requirements

Active cavity control is enabled if the following requirements are both met:

- APC (Constant Power) mode
- Power > 0.8 mW


To disable/enable active cavity control

To disable active cavity control, user simply press the **2nd** + **APC** key combination (the APC LED is blinking). To toggle between *active control* On/Off, press the **2nd** + **APC** command keys again.

System status is summarized in the following table:

APC LED	Constant Power Mode	Active Cavity Control
On	Yes	Yes
Blinking	Yes	No ⁹
Off	No	No

Note: If you press the APC key and the APC LED went off, this indicates that the system is now in constant current mode. When entering constant current mode, the active cavity control becomes disabled, and the system restores the operating parameters that were in effect at the time TUNICS last left the constant current mode.

	IMPORTANT	Active Cavity Control & External Modulation Prior to implementing external modulation via “LF Mod”, “HF Mod”, and “FSC” input connectors, user must first disable the system’s active cavity control by pressing the 2nd + APC keys.
---	------------------	---

⁹

If the APC LED starts blinking while the TUNICS-**Plus** is in APC mode with active cavity control on, this may indicate that the laser power output (P) is set lower than the 0.8 mW threshold.

Active Wavelength Monitoring

When associated with an external optical wavelength meter, the TUNICS-**Plus** tunable laser can implement active wavelength monitoring to achieve high-precision wavelength tuning (± 2 pm, typically; or within the wavemeter's accuracy).

Instrument Setup and System Requirements

The instrument setup and system requirements are outlined below:

1. The measurement setup requires optical connection of the TUNICS-**Plus** output to the wavemeter via the low output port of an optical coupler (the remainder of the signal available on the coupler's high-output is sent to the component under test). (See Figure 5 - 1.)

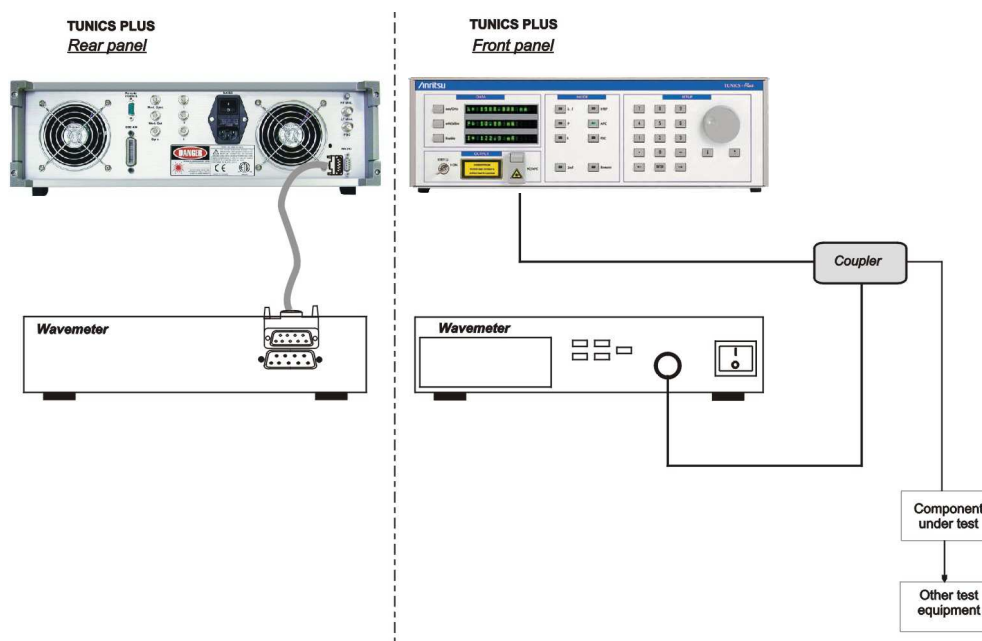


Figure 5 - 1: Instrument Setup for Wavelength Monitoring

2. Active wavelength monitoring requires a communication link between the TUNICS-**Plus** tunable laser source and the wavemeter via the RS-232 C interface. Make sure the RS-232 C cable connection meets the following pin assignment on the TUNICS-**Plus** SubD-9 connector:
 - Pin 2: RX (Data received by TUNICS-**Plus** on pin 2)
 - Pin 3: TX (Messages sent out by TUNICS-**Plus** via pin 3)
 - Pin 5: GND (The ground is through pin 5)

The wavemeter you are using must meet the following hardware and configuration requirements:

Serial port	
Command set	: Based on the SCPI (Standard Commands for Programming Instruments) syntax.
Serial interface	
configuration	: RS-232 → ON
	Baud rate → 9600
	Handshaking → None

Depending on the wavemeter you are using, configuration can be achieved either via soft-panel controls or by hardware setting via a backpanel dip-switch. For details on how to configure your wavemeter, please refer to the user's manual provided by the manufacturer.

Note: as an example, with a Burleigh-EXFO wavemeter, it is necessary to use a DB9 male-male, straight-through serial cable, to connect the TUNICS-**Plus** RS-232 C port to the wavemeter serial port.

(Please contact us for a list of compatible wavemeters.)

Wavelength Monitoring Scheme

The wavelength monitoring scheme consists in a loopback to the tunable laser of the wavemeter actual reading, and self-alignment of the tunable laser by iteration process to achieve the user-defined target wavelength within the wavemeter's accuracy. This is outlined in the following diagram:

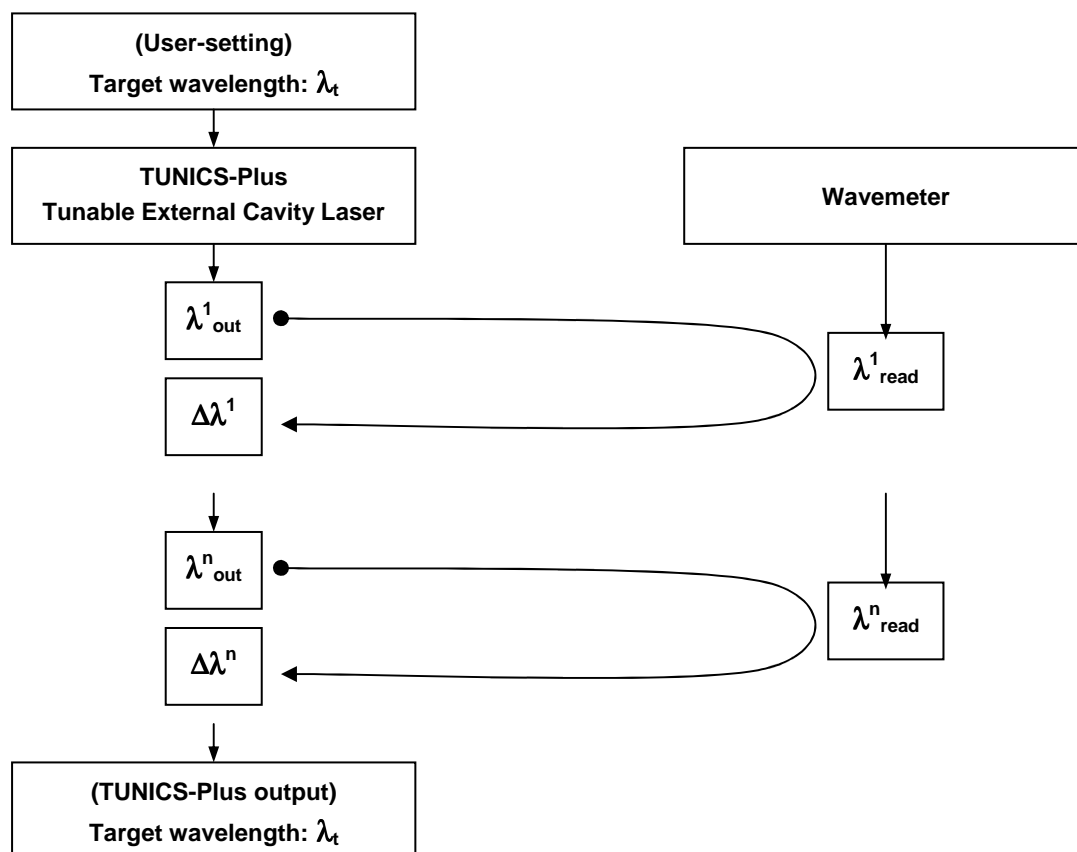


Figure 5 - 2: Wavelength Monitoring Scheme

Depending on the performance of the wavemeter you are using, the average iteration time to reach the target wavelength may range from 1.5 s to 2 s.

How to Perform Wavelength Monitoring

To perform wavelength monitoring of the TUNICS-**Plus** using a wavemeter, follow these steps:

Step	Action
1.	Configure the wavemeter remote interface with the following parameters: <ul style="list-style-type: none"> • Remote: [On] • RS-232 : [On] • Baud: [9600] • Handshaking: [None]
2.	Connect the TUNICS- Plus DB9 serial port to the wavemeter serial port.

3.	As the tunable laser acts as the controller device for the wavemeter, press the 2nd + Remote function keys at the TUNICS- Plus front-panel and select: • RS-232 IN : [OFF] ¹⁰
4.	Start the Feedback function by pressing the 2nd + dBm/mW keys.
5.	Enter the wavelength setting at the TUNICS- Plus front-panel.
<u>Result:</u>	The TUNICS- Plus executes the wavelength monitoring scheme (see Figure 5 - 2) so the wavelength of the output signal matches closely the wavelength setting. The message “ Moving... ” is displayed at the front-panel indicator till the output wavelength has reached the target value and has stabilized within wavemeter accuracy.
<u>Note:</u>	To compensate for possible disturbance such as accidental vibrations, shocks, temperature drifts, and possible malfunction of the wavemeter, TUNICS- Plus performs the wavelength monitoring scheme with a 4 seconds periodicity.
<u>Display:</u>	The TUNICS- Plus front-panel display shows the tunable laser operating wavelength, alternating with the “λ FEEDBACK:ON” message to indicate that the “active wavelength monitoring” is currently on. Immediately after you turn off the “active wavelength monitoring”, the message “λ FEEDBACK:OFF” shows once then goes off. If the “active wavelength monitoring” does not work properly, the “λ FEEDBACK:ERR” message shows in alternance with the tunable laser operating wavelength.

¹⁰

While in wavelength monitoring configuration, the TUNICS-**Plus** tunable laser can only be remote operated via the IEEE-488.1 interface. If you wish to restore remote operation of the TUNICS-**Plus** via the RS-232 port, you must press the **2nd** + **Remote** keys at the front-panel and select, RS-232 IN : [ON].

Auxiliary Inputs and Outputs

Figure 5 - 3 shows the auxiliary input and output connectors on the rear panel of TUNICS-Plus.



Figure 5 - 3: Auxiliary Input/Output Connectors on Rear Panel of TUNICS-Plus

In addition to the two remote control connectors (RS-232 C and IEEE-488.1), eight BNC connectors (λ , I, P, Mod Sync, Mod Out, Sync., LF Mod., and FSC), one SMA connector (HF Mod.), and one remote interlock connector are also available on the rear panel of the system.

External Modulation

	IMPORTANT	Active Cavity Control & External Modulation Prior to implementing external modulation via “LF Mod”, “HF Mod”, and “FSC” input connectors, user must first disable the system’s active cavity control by pressing the 2nd + APC keys.
--	------------------	---

Intensity Modulation - Low Frequency

It is possible to implement low frequency intensity modulation of the output beam from DC to 8MHz in the Direct Current mode and from 30 kHz to 8 MHz in the APC mode.

Modulation of the diode current is obtained by applying a signal anywhere between **0 and 5 volts** on the BNC connector labeled “LF Mod.”

LOW FREQUENCY MODULATION SETUP

To implement low frequency modulation, user can connect an external modulation source to the “LF Mod.” BNC connector, as shown in Figure 5 - 4.

Note: The output impedance of your waveform generator must match the 50 Ω input impedance of the TUNICS “LF Mod.” BNC connector. Also, make sure to use a 50 Ω BNC cord to connect the TUNICS and the waveform generator.

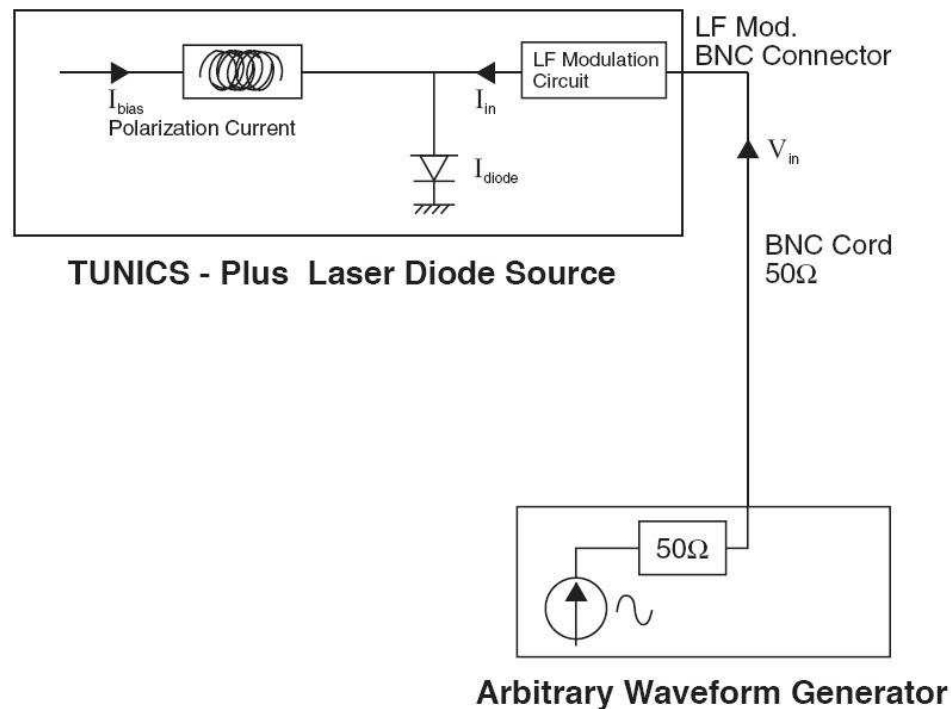


Figure 5 - 4: Low frequency modulation setup

USING LOW FREQUENCY MODULATION

User first set the operating current I_{bias} of the TUNICS in constant current mode. The relationship between the diode output current I_{diode} and the input voltage V_{in} is then given by:

$$\text{Formulae: } I_{diode} = I_{bias} - \frac{V_{in}}{R_{internal\Omega}}$$

Where:

- V_{in} is the input voltage applied to LF Mod. BNC input connector (V_{in} must always be positive and less than 5 V)
- I_{diode} is the total diode current
- I_{bias} is the polarization current of your TUNICS
- $R_{internal\Omega}$ is the 50Ω resistor of the LF modulation circuit

With a maximum allowed voltage input of 5 Volts, the maximum input current that user can add to the TUNICS laser diode is 100 mA.

RESULTING LF-MODULATED DIODE CURRENT

The laser diode of your TUNICS is characterized by its saturation current (I_{max}) and the cutoff current (I_{cutoff}). To perform a diode intensity modulation under optimum conditions, you must ensure at all times that the resulting modulated diode current obeys the following condition: $I_{cutoff} \leq I_{diode} \leq I_{max}$.

Operation

Refer to Figure 5 - 5 below for a behavioral representation the resulting modulated diode current.

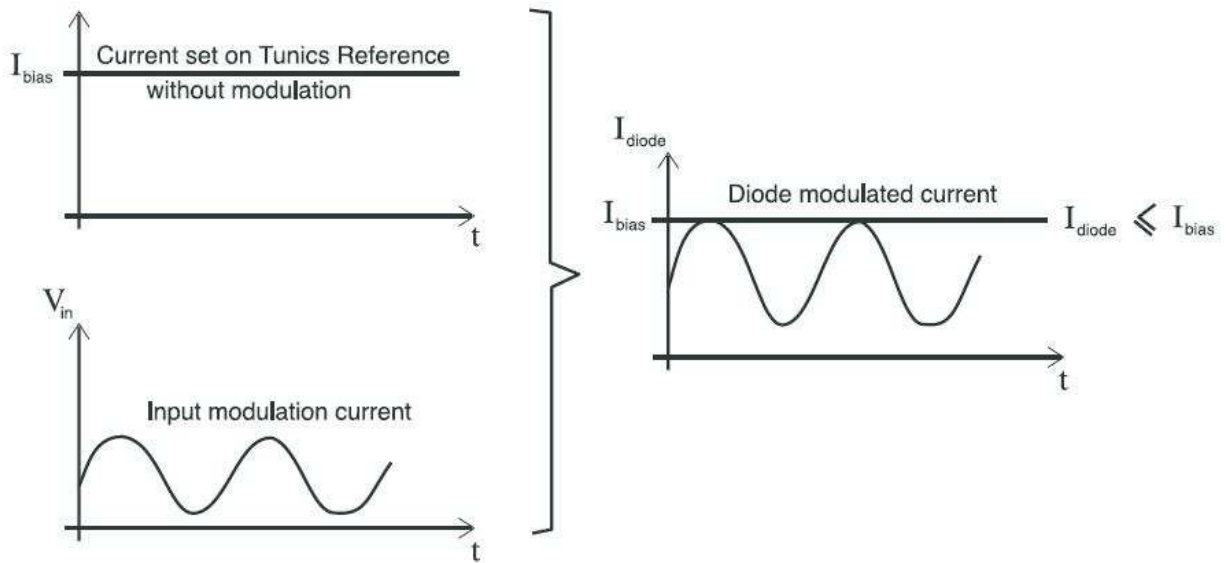


Figure 5 - 5: Resulting LF-modulated diode current

Note: In the APC mode, the output power is controlled in order to remain constant, which is in principle not compatible with a modulation of the output power. However, if the frequency of the modulation signal is well above the intrinsic frequency ($\cong 1$ kHz) of the APC loop, the average power is monitored, and remains constant. Modulation is possible in the APC mode at frequencies above 30 kHz. If a lower frequency is used, interaction between modulation signal and APC circuit will cause unpredictable modifications of the output power waveform.

Intensity Modulation - High Frequency

The “HF Mod” input allows user to increase the total diode current by applying voltage to the “HF Mod” SMA connector. A valid high frequency modulation signal applied to the SMA port can range from 30 KHz up to 1 GHz.

HIGH FREQUENCY MODULATION SETUP

To implement high frequency modulation, user can connect an external high frequency modulation source to the “HF Mod.” SMA connector, as shown in Figure 5 - 6.

Make sure to use a 50Ω coaxial cord to connect your HF generator to the 50Ω input SMA connector.

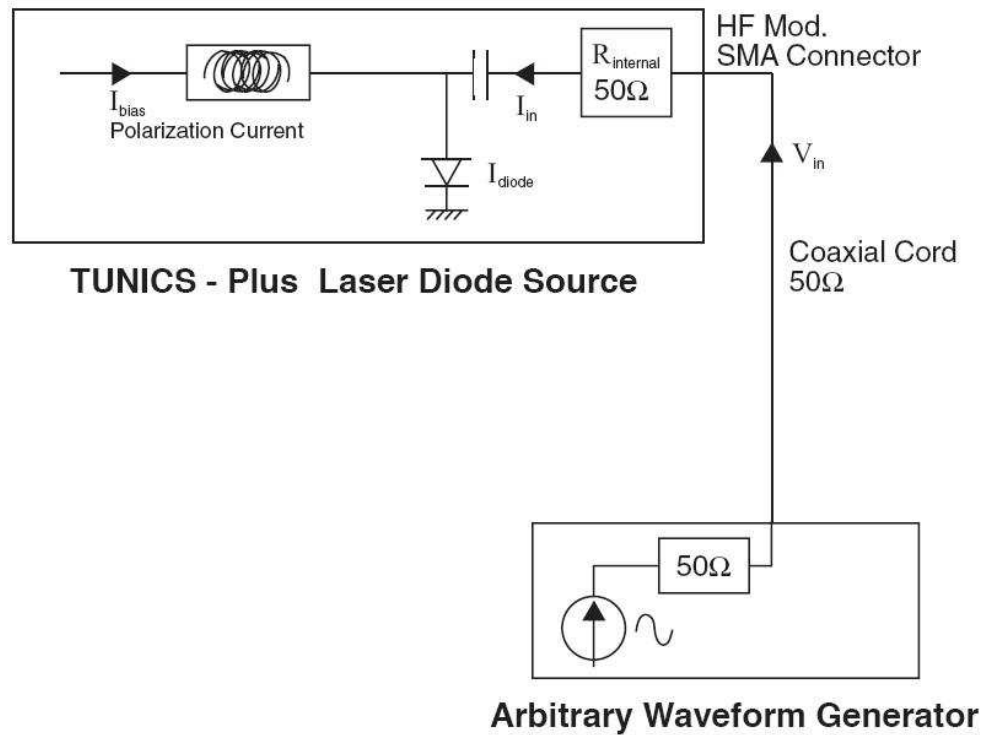


Figure 5 - 6: High frequency modulation setup

USING HIGH FREQUENCY MODULATION

The relationship between the voltage at the “HF Mod” input and the current added to the diode is given by:

$$I_{in} = V_{in}/R_{internal}$$

where: I_{in} is the added current

V_{in} is the voltage applied to the “HF Mod” input

$R_{internal}$ is the 50Ω internal resistor of your TUNICS

By modulating this extra current, it is possible to modulate the optical power output within a frequency range from 30 kHz to 1 GHz. The low-frequency cutoff is 30 kHz. The modulation bandwidth is greater than 1 GHz, and this “HF Mod” input can be used for mode-locked operation as well.

	WARNING	<p>Since the “HF Mod” input is directly connected to the laser diode chip, user must ensure that the input voltage on the “HF Mod” connector is such that the following conditions are always met:</p> <ol style="list-style-type: none"> 1. The diode laser is forward biased. 2. The current flowing through the diode is lower than the maximum admissible current (I_{max}) indicated in the system acceptance test report.
--	----------------	--

The best way to fulfill the above conditions is to measure the system bias current (I_b) through the voltage on the current monitor output, and to compute the total current using the following formula:

$$I_{diode} = I_{bias} + V_{in}/R_{internal}$$

Where:

- I_{diode} is the total current through the diode
- I_{bias} is the bias current
- V_{in} is the voltage applied to the “HF Mod” input
- $R_{internal}$ is the 50Ω TUNICS internal resistor

RESULTING HF-MODULATED DIODE CURRENT

User must monitor V_{in} to ensure that the total current remains positive and lower than the maximum admissible current:, as shown in Figure 5 - 7 below:

$$0 < I_{diode} < I_{max}$$

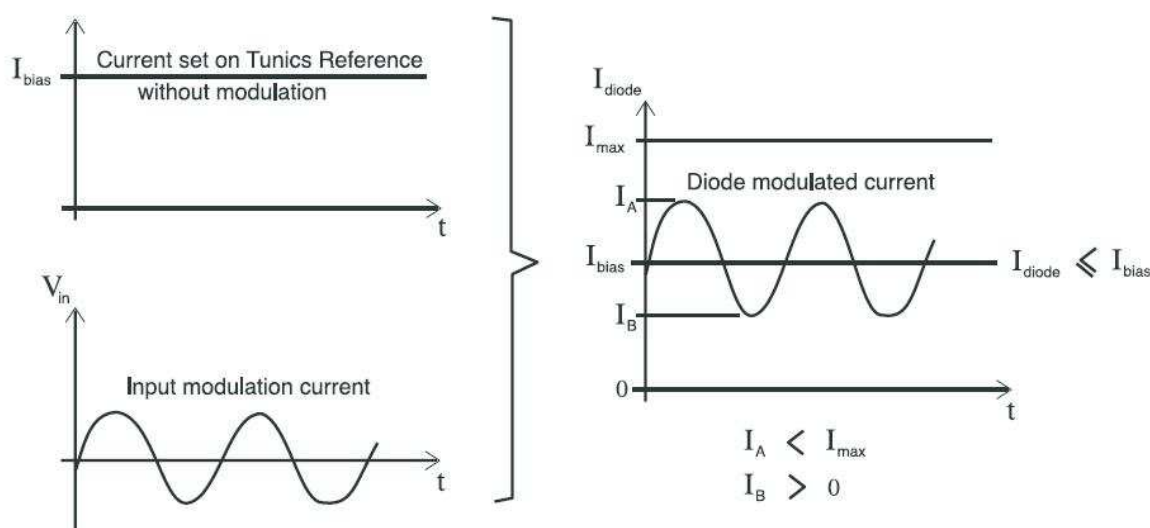


Figure 5 - 7: Resulting HF-modulated diode current

If the voltage applied to the “HF Mod” input is an AC signal with no DC component, and whose peak-to-peak amplitude is V_{pp} , the above condition becomes:

$$V_{pp} < 2 \times R_{internal} (I_{max} - I_{bias})$$

$$V_{pp} < 2 \times R_{internal} (I_{bias})$$

Therefore, V_{pp} must be lower than the lowest of these two values.

Note: TUNICS constantly monitors the average optical output power and regulates the system bias current to keep the system output power constant. If a DC or low frequency (below 30 kHz) voltage is applied to the “HF Mod” input, there is no modulation on the output power, as the system automatically compensates by adjusting the bias current to cancel this modulation. However, if the frequency of the signal applied to the “HF Mod” input is well above the power control loop cutoff frequency (30 kHz), the output signal modulation will match that of the input signal, around an average power determined by the settings of the system. For that reason, user must implement a modulation frequency above 30 kHz.

	CAUTION	The “HF Mod” SMA input is directly connected to the laser diode chip via a 50Ω resistor. It is essential that the user insure that at all times the diode is forward biased to avoid the destruction of the diode.
--	----------------	--

Internal Modulation

A 7.2 kHz TTL modulation signal is provided at the rear panel of the electronics controller. The corresponding BNC connector is labeled "Mod Out". This TTL modulation signal can be used to digitally modulate the TUNICS by connecting back to the "LF Mod" BNC connector.

Synchronization Signal

A 7.2 kHz TTL synchronization signal is also provided on the BNC connector labeled "Mod Sync".

Wavelength Modulation

WAVELENGTH FINE SCANNING

TUNICS-**Plus** has a wavelength fine tuning capability which can be used to adjust the wavelength very precisely, or to serve the optical frequency. This function is available via the FSC connector on the rear panel. When voltage is applied to the FSC connector, the optical frequency shifts 300 MHz for every volt applied. The maximum input voltage allowed on FSC is ± 10 volts. FSC input also allows wavelength modulation, by applying a modulation signal whose frequency is lower than 10 kHz.

COHERENCE CONTROL

When a "noise" signal is directly fed into the "FSC" input via the **2nd** + **FSC** function to modulate the optical frequency, this induces an apparent broadening of the laser linewidth leading to a reduction of the coherence length. This can be useful to eliminate parasitic interferometers in the user setup. This modulation is limited to low frequency (DC to 10 kHz typ.). When this function is active, the LED on the FSC key is blinking.

Monitor Outputs

Two output connectors located on the rear panel enable continuous monitoring of actual values of current (**I** output), power (**P** output), wavelength (**λ** output), and detection of motor displacement (**Sync.** output). These monitor outputs use BNC connectors. Signal levels are as follows:

Monitor Output	Description	Reading
I output	Current monitoring	20 mV/mA (typ)
P output	Power monitoring	0.3 mV/ μ W (typ)
λ output	Wavelength monitoring	See table below
Sync.	Power monitoring	5 V: motor moving 0 V: motor idle

λ output is summarized in the following table:

TUNICS Model	Plus O	Plus E	Plus S
λ for V=0 (nm)	1230	1320	1400
λ output (mV/nm)	33	33	33

TUNICS Model	Plus S WB	Plus CL	Plus CL WB
λ for V=0 (nm)	1360	1490	1460
λ output (mV/nm)	25	33	25

Remote Interlock

TUNICS-**Plus** is a Class IIb product (21 CFR Subchapter J) equipped with a Remote Interlock connector. The user can attach a remote interlock switch to the Remote Interlock SUBD-9 connector located on the rear panel.

When connecting the remote interlock switch, user must provide a connector that meets the following pin assignment:

Remote Interlock Connector Pin Assignment

Pins 1-2 Shorted
Pins 6-7 Connect to external switch

Note: User must ensure the external switch is electrically isolated from other circuits, including earth ground.

Opening the interlock switch disables the laser output. The laser output may be turned on again by simply closing the interlock switch, and pressing the **Enable** key.

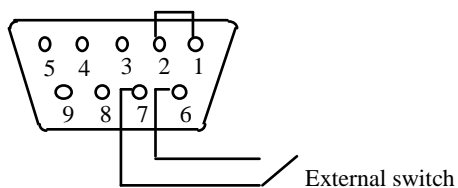


Figure 5 - 8: Remote Interlock SUBD-9 connector

REMOTE CONTROL

Introduction

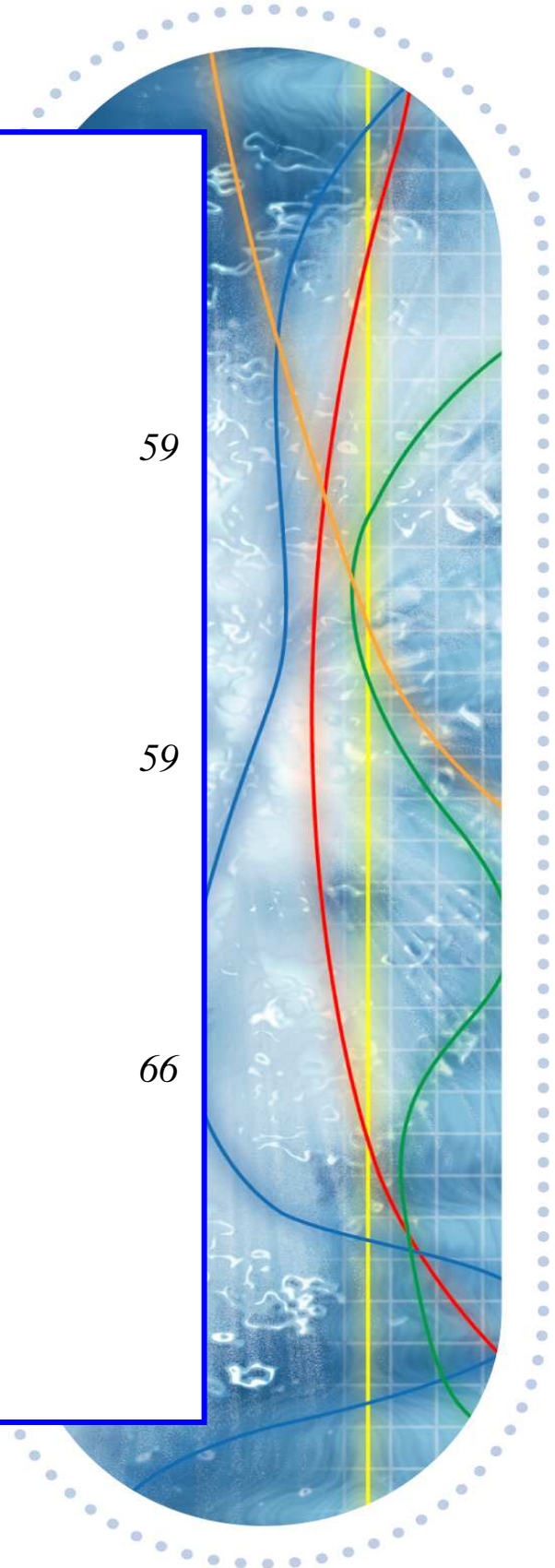
59

RS-232 C Remote Control

59

IEEE-488.1 Remote Control

66



Remote Control

Introduction

TUNICS-**Plus** features two types of remote control interfaces, RS-232 C and IEEE-488.1 connectors on the rear panel. When TUNICS-**Plus** receives data from either the RS-232 C or the IEEE-488.1 interfaces, it automatically enters the remote control mode and the LED on the Remote key is lit to indicate that remote control is active. When TUNICS-**Plus** is operating in remote mode, all function keys on the front panel are disabled at the exception of the Remote key, which allows user to switch between remote and manual modes. To exit the Remote mode without breaking the cable connection to the remote terminal or computer, simply press the **Remote** key.

RS-232 C Remote Control

To operate in RS-232 C remote mode, connect a VT-100 compatible terminal or a PC equipped with a terminal emulation program to the RS-232 C port on the rear panel of TUNICS-**Plus**.

Note: you first need to activate the RS-232 remote control mode, as explained in section Front & Rear Panels, in the Mode Area, page 28.

When the Remote key LED is lit, all front panel manual commands are rerouted to the remote terminal or computer via the RS-232 C interface. In addition, user interface specific commands are also available in remote mode.

Physical Interface

The cable link is a crossover, three-wire RS-232 C type that connects the TUNICS-**Plus** SubD-9 connector to a computer serial port. Make sure the RS-232 C cable connection meets the following pin assignment on the TUNICS-**Plus** SubD-9 connector:

- Pin 2: RX (Data received by TUNICS-**Plus** on pin 2)
- Pin 3: TX (Messages sent out by TUNICS-**Plus** via pin 3)
- Pin 5: GND (The ground is through pin 5)

Note that no other connector-pin is used.

Make sure to configure the computer's terminal emulation program with the following parameters:

Data transmission rate	9600 bauds
Number of bits	8
Parity	None
Stop bits	1

The following table shows how to connect the TUNICS-**Plus** system (DCE) to your PC (DTE):

Table 6 - 1: Crossover, 3-wire RS-232 serial connection to TUNICS- **Plus**

Digital Communications Equipment (DCE)			Digital Terminal Equipment (DTE)		
1	o			o	1
(RXD) 2	o			o	2
(TXD) 3	o			o	3
4	o			o	4
(GND) 5	o			o	5
6	o			o	6
7	o			o	7
8	o			o	8
9	o			o	9

Message Format

All commands sent by the computer to the TUNICS-**Plus** consist of a string of ASCII alphanumeric characters followed by a “carriage return” character (CR, or ASCII code 13). When TUNICS-**Plus** receives such an instruction, the character string is decoded and the relevant procedure is carried out.

After completion, TUNICS-**Plus** returns a message to inform the computer that the requested procedure has been performed or that an error event has occurred. This response message is always ended with the same end-of-message group: “↵” (carriage return), “>” (greater than), and “ ” (space character).

These three characters indicate that the TUNICS-**Plus** system is ready to receive a new instruction.

Table 6 - 2: Example of dialog between the computer and TUNICS-Plus** (RS-232 C)**

Computer/Terminal → TUNICS- Plus		TUNICS- Plus → Computer/Terminal	
Command String	Instruction	Response Message	Description
APCON	Set Constant Power (APC) On.	OK↵ >	Instruction performed correctly. TUNICS is awaiting another instruction.
P = 30↵	Change power output value to 30 mW.	Value error↵ >	Power output is over the maximum allowed. TUNICS is awaiting another instruction.
L?↵	What is the present wavelength value?	L = 1523.32↵ >	Present wavelength value. TUNICS is awaiting a new instruction.

The message syntax described above eases the management of the TUNICS-**Plus**/Computer interface since the computer must always await the acknowledgement before sending a new instruction. In addition, the response from TUNICS-**Plus** is

always ended with the three characters “↵”, “>” and “space”, enabling user to easily spot the end of a message.

Note 1: TUNICS-**Plus** user interface commands are not case sensitive.

Note 2: Multiple commands may be sent in single line and should have instructions separated by a semicolon “;”, and a command string must always end with a carriage return character “↵”. The various instructions are processed by TUNICS-**Plus** as soon as the “↵” character is entered and executed in the order corresponding to the command string sequence (make sure sequential instructions are logically placed within a command string). However, responses generated by multiple command strings are ended with the carriage return character “↵”. Only the last response is ended with the end-of-message group “↵”, “>” and “space”.

Note 3: The TUNICS-**Plus** input buffer is 255 characters long. Therefore, if you enter a command string more than 255 characters long, or submit a new command before previous strings have been processed, the buffer will be cleared, all received instructions lost, and the message “command error” will be sent.

Note 4: “White space”¹¹ characters are allowed in some places inside the command strings sent to TUNICS-**Plus**. They can be placed at the beginning of an instruction, after an instruction, as well as before, after, or in place of the “=” sign. They cannot be inserted within an instruction mnemonic, within a numeric value, or between a mnemonic and the “?” mark.

Response Messages

All successful instructions resulting in a parameter change or a change of operating mode produce an “OK” statement ended with “↵”, “>”, or “space”.

If a command string is not recognized by TUNICS-**Plus**, the message returned will be “**COMMAND ERROR**”. If a command requesting a new parameter setting is correctly formulated, but the specified value is outside the allowed limits for that parameter, the setting will remain unchanged and the response is “**VALUE ERROR**”.

Format of Numerics

Some command messages include a value. The format rules for these values are as follows:

- No spaces (see *Note 4* above) are allowed within values; leading zero (“0”) characters are allowed at the beginning of a value.
- Digits after the decimal point may be sent or omitted.
- Digits after the decimal point are sent after the characters “.” or “,”.
- Values are ended with the “CR” character.
- Do not specify any units to parameter values. All values take the default units, factory-set for TUNICS-**Plus**.

¹¹

On the RS232 C interface, all characters with ASCII codes lower than or equal to 32, except the carriage return character (ASCII code 13) are considered as white space characters.

Table 6 - 3 provides some correct and incorrect examples of the message formats for numerics.

Table 6 - 3: Examples of message formats with numerics (RS-232 C)

Entry	Correct/Incorrect	Meaning
I=5	Correct	Set diode current to 5.0 mA
P=0.22	Correct	Set power to 0.22 mW.
L=1530.2	Correct	Set wavelength to 1530.200 nm.
P=01	Correct	Set power to 1 mW.
P= 0.65	Correct (spaces before the numeric are allowed.)	Set output power to 0.65 mW.
P=0.65 mW	Incorrect (unit abbreviations are not allowed.)	Command error
L=1 520.31	Incorrect (spaces within the numeric between “1” and “5” are not allowed.)	Command error

Changing Parameter Settings

When “nn.n” occurs below, it stands for the transmitted numerical value. All commands must end with the “CR” character.

Responses sent by TUNICS-Plus may be as follows:

OK	Command successfully executed
VALUE ERROR	Value outside valid range
COMMAND ERROR	Syntax error

Table 6 - 4 lists the commands for each parameter setting of the tunable laser.

Table 6 - 4: Command parameters (RS-232 C)

Command	Action
I=nnn.n	Sets the laser current level (in mA) and switches to the constant-current mode.
P=nn.nn ¹²	Sets the optical power (in mW) and switches to the constant-power mode.
P=±nn.nn ¹²	Sets the optical power (in dBm) and switches to the constant-power mode.
L=nnnn.nnn	Sets the emission wavelength (in nm). The laser is tuned directly to the desired wavelength without step-mode scanning or swept-mode hop-free scanning. An “OK” statement is received when the emission wavelength has stabilized.
F=nnnnnn.n	Sets the optical frequency (in GHz). The laser is tuned directly to the desired optical frequency without step-mode scanning or swept-mode hop-free scanning.. An “OK” statement is received when the emission wavelength has stabilized.
FSCL=nn.n	Switches to the Fine Scanning mode (canceled by any other command). “nn.n” is the wavelength change in picometers (pm).
FSCF=n.nn	Switches to the Fine Scanning mode (canceled by any other command). “n.nn” is the optical frequency change in GHz.

¹²

The units (dB or mW) should be first selected with the appropriate DBM or MW command.

Table 6 - 4: Command parameters (RS-232 C) (cont'd)

B_SUPPR=n	0=Eliminates the effect of the tuning micrometer screw backlash. 1=Restores the backlash suppression procedure.
MOTOR_SPEED=nnn	Sets the sweep speed (in nm/s) used in swept-mode, hop-free scanning operation. Valid speed settings can range from 1 to 100 nm/s. Operational sweep speeds are: [1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 17, 18, 20, 22, 25, 29, 33, 40, 50, 67, 100] nm/s. Therefore, the speed setting you provided is automatically rounded to the nearest operational sweep speed. For example, if you provide "MOTOR_SPEED=060", the actual sweep speed implemented by the system is 67 nm/s. To check the current operational sweep speed, use the "MOTOR_SPEED?" query.
PCAL1 =nn.nn	Sets the calibration power (@ first wavelength).
PCAL2 =nn.nn	Sets the calibration power (@ second wavelength).

Changing the Operating Mode

Table 6 - 5 lists the commands that change the TUNICS-**Plus** operating mode.

Table 6 - 5: Operating mode commands (RS-232 C)

Command	Action
APCON	Switches system operation to the Automatic Power Control (APC) also known as "constant-power" mode.
APCOFF	Cancels the APC and switches system operation to the "constant-current" mode. Note: Using the APCON and APCOFF commands is the same as pressing the APC key in manual operating mode to toggle between constant power and constant current mode.
ACTCTRLON	Enables "active cavity control". (For more details on "active cavity control", please refer to page 45 in this manual.)
ACTCTRLOFF	Disables "active cavity control".
AUTO_CAL	Launches the system auto-calibration sequence. (Same as pressing the 2nd + I keys in manual mode.) Note 1: Using either the AUTO_CAL or INIT commands forces the system into performing an auto-calibration. Note 2: If the auto-calibration has failed, the message "referencing error" is displayed for a few seconds, and the TUNICS resumes normal operating status. Unlike the 2nd + I command that stalls the system when a referencing error occurs, there is no need to press Enter to restore the system.
DBM	Sets dBm as the unit for optical power values.
MW	Sets mW as the unit for optical power values.
INIT	Starts the initialization of the optical head (this takes at least 10 s). Note: The INIT command is available in RS-232 C mode only. . It generates the answer "Software Vn.nn" which stands for the software version. (Using this command in IEEE-488 mode would cancel all the pending commands stored in the GPIB buffer.)
ENABLE	Enables the optical output.
DISABLE	Disables the optical output.

Table 6 - 5: Operating mode commands (RS-232 C)

L_FEEDBACK=1	Enables the “active wavelength monitoring” while using an external wavemeter (see also page 46).
L_FEEDBACK=0	Disables the “active wavelength monitoring”.
CTRLON	Sets the linewidth coherence control ON.
CTRLOFF	Sets the linewidth coherence control OFF.

Note: Wavelength scans can only be performed manually, via the command keys at the TUNICS-**Plus** front panel.

Reading System Parameters

Table 6 - 6 lists the queries for system parameters of the TUNICS-**Plus** tunable laser.

User can monitor system parameters such as laser output power, current level, and emission wavelength at any time during operations by typing appropriate parameter queries. A valid parameter value query consists of the command name immediately followed by the “?” mark.

Table 6 - 6: System parameters (RS-232 C)

Query Statement	Action	System Response
I?	Returns the present diode-current level.	“I=nnn.n” (in mA) or “disabled” if the Enable mode is not active.
P?	Returns the present laser output power level.	“P=nn.nn” (in mW) or “P=±nn.nn” (in dBm) or “disabled” if the Enable mode is not active.
L?	Returns the present emission wavelength value.	“L=nnnn.nnn” (in nm).
F?	Returns the present optical frequency value.	“F=nnnnnn.n” (in GHz).
MOTOR_SPEED?	Returns the operational sweeping speed.	“MOTOR_SPEED=nnn” (in nm/s) (1 ≤ nnn ≤ 100)
LIMIT?	Returns the status of the current limitation feature.	“Yes” if the current is limited and “No” otherwise.
B_SUPPR?	Returns the status of the backlash suppression.	0 = the backlash suppression procedure is not in use. 1 = the backlash suppression procedure is in use.
PCAL1?	Returns the calibration power (@ first wavelength).	“PCAL1=nn.nn” (in mW).
PCAL2?	Returns the calibration power (@ second wavelength).	“PCAL2=nn.nn” (in mW).
ENABLE?	Returns the state of the laser output.	“ENABLE” or “DISABLE”
L_FEEDBACK?	Indicates whether the “active wavelength monitoring” function is in use.	“1”: in use “0”: turned off “ERROR”: possible malfunction (See <i>Appendix A: Troubleshooting</i> for details.)
*IDN?	Returns system ID.	

Other System Commands

Table 6 - 7 lists the other system commands available on the TUNICS-**Plus** tunable laser.

Table 6 - 7: Other system commands (RS-232 C)

Command	Action
ECHON	Activates "echo" mode. When this mode is active, TUNICS sends an echo of each character received back through the serial cable. Since some terminals and terminal emulation programs do not feature local echo, this function can be useful for visual monitoring of the characters keyed in at the terminal or for having a secure link.
ECHOFF	Cancels "echo" mode. Note: The default system setting is ECHOFF. In addition, ECHON mode is restored to ECHOFF when user presses the Remote key to return to manual-mode operation.
LOCAL	Switches back to local mode.

IEEE-488.1 Remote Control

GPIB Overview

The "GPIB" (General Purpose Interface Bus) communications link of TUNICS-Plus is compatible with the IEEE-488.1 Standard. The purpose of the normalized IEEE-488 bus is to allow communications between instruments and computers through a specialized parallel interface. The data transfer rate is as high as 1 Mbytes/s, and the information transfer is highly secured. The GPIB link is often used to automate measurement processes.

The standard configuration is composed of a controller device (which is usually a computer equipped with a GPIB interface board and corresponding IEEE-488 terminal emulation software) linked to various instruments. The controller manages the flow of information to, from, and between devices. One "talker" device and one or more "listener" devices may be present at any time on the GPIB link. Data sent on the bus is encoded in ASCII strings.

Each instrument on the bus is identified by its own GPIB address. The GPIB address is an integer in the range of 1 to 30. For each instrument, different control registers can be used to determine the instrument operation state. The default TUNICS-Plus GPIB address is 10 but it can be changed (see *Changing the GPIB Address* on page 71).

GPB Functions Available on TUNICS-Plus

TUNICS-Plus supports the following GPIB functions, as listed in the table below:

Table 6 - 8: Supported GPIB functions

Mnemonic	Meaning	Implementation
SH1	Source handshake	Complete
AH1	Acceptor handshake	Complete
T5	Talker	Complete
L3	Listener	Complete
SR1	Service request	Complete
RL1	Remote / Local	Complete
PP0	Parallel poll	No
DC1	Device clear	Complete
DT0	Device trigger	No
C0	Controller	No

Checking System Status with the Status Word Variable

Definition of the Status Word

The status word is an 8-bit value that reflects the status of the TUNICS-**Plus**. It contains a number of binary indicators, which can be used by the computer/terminal for optimal synchronization with the TUNICS-**Plus** unit. They indicate to the computer/terminal the nature of the current operations as well as the errors encountered. The only way to read the status word is to perform an operation called serial polling.

Table 6 - 9: Bit layout in status word (GPIB)

Bit #	Mnemonic	Meaning
7 (MSB)	(Bit not used)	Bit 7 is not used.
6	SRQ (Service ReQuest)	Set when a service request has been sent by the instrument. This bit is cleared when the bus controller reads the status register.
5	ESB	(Reserved for future use.)
4	MAV (Message AAvailable)	Set when a message is placed in the output buffer. This bit remains set for as long as the output buffer is not polled.
3	LIM (Current LIMitation)	Activated when TUNICS- Plus is current-limited while in APC mode. Deactivated in all other cases.
2	ERRV (ERror in Received parameter Value)	Indicates that the value entered as a parameter could not be read, exhibits an invalid format, or is outside the valid range (e.g.: I=160). (For example, a parameter value becomes invalid when the parameter's unit is provided.) (The indicator is disabled when a correct instruction is received.)
1	ERRC (ERror in Received Command)	Indicates that an invalid command has been received. (Disabled when a valid command is received.)
0	OPC (OPeration Complete)	Indicates that the execution of the last command is complete. TUNICS- Plus is ready to receive new instructions. Activated when no task is currently performed by TUNICS- Plus . Deactivated when handling commands (and when shifting the drive to modify the wavelength emission). Particularly, during a wavelength scan operation, the OPC bit-value is 0 each time the motor is moving and 1 during the pause between two successive steps.

Note: when a condition is in effect, the corresponding binary indicator takes a bit-value of (1); otherwise, the default bit-value is (0).

Use of the Service Request (SRQ)

The SRQ (Service Request) line is part of the definition of the IEEE-488.1 standard. The SRQ line can be triggered by any device on the bus. It alerts the central controller that a particular device requests to perform a specific operation.

When the controller detects that the SRQ line is active, it can probe all of the devices present on the bus to determine which device has initiated the service request and for what reason. This operation by which the controller reads the status word of each device is called “serial polling”. The controller then monitors the bit-value of each indicator within the status word (see below).

The IEEE-488.1 standard specifies that when a device sends a service request, bit-value # 6 is set to 1. The other bits of the status word can reflect the state of different logical indicators of the system.

To help synchronize operations between the TUNICS-**Plus** and the computer/terminal, the latter can instruct the TUNICS-**Plus** to send a service request whenever certain bit-values are set to 1, as defined by the user. This can be achieved by using the instruction *SRE. Therefore, user does not need to repeatedly prompt for the status word since the TUNICS-**Plus** automatically notifies the computer/terminal when a particular event occurs.

User can configure the *SRE instruction to perform a particular service request. The transmitted parameter associated with the instruction *SRE is an integer that takes a value between 0 and 255. The *SRE can also be presented as an 8-bit binary value, each bit-value is either 0 or 1 depending on whether the corresponding status indicator is off or on.

For example, if the instruction *SRE=16 is received by the TUNICS-**Plus**, this corresponds to the binary value 00010000, as shown in Table 5-12. Since in this example the bit indicator #4 is set to 1 (Bit #4 = MAV = Message Available), the SRQ line will be automatically activated each time a message becomes available.

The following table provides examples of SRE values when individual indicator bit-values are set to 1.

Table 6 - 10: Examples of single-indicator SRE values

Indicator	SRE decimal value	Bit position #	Binary value
MAV	16	indicator #4	00010000
LIM	8	indicator #3	00001000
ERRV	4	indicator #2	00000100
ERRC	2	indicator #1	00000010
OPC	1	indicator #0	00000001

To calculate the parameter value to send with the instruction *SRE, simply add together the decimal values found in Table 5-12 of each individual indicators you would like to combine. For example, to obtain a service request each time an error is made, you must send the instruction “*SRE=6”, obtained by adding the corresponding values of the two error indicators, which are ERRV(4) and ERRC(2).

Once the TUNICS-**Plus** has activated a service request, the status word is no longer automatically updated, until the controller performs a new serial poll on the GPIB link.

As the status word indicators remain unchanged, this allows the status word to provide reliable system-status information at the time of service request.

IMPORTANT: Since the processing of an instruction begins with the deactivation of the OPC indicator, no new instruction can be processed until the computer has performed the serial poll. Therefore, user should perform a serial poll as soon as a service request is received by the computer. Once the status word is read, the service request is deactivated and the TUNICS-**Plus** can resume normal operation.

Syntax for Parameter Commands

Syntax and Message Conventions

Each data is sent over the bus as an ASCII character string. The following syntax applies to message send and receive:

- Each message string sent by the computer/terminal to the TUNICS-**Plus** is ended with the “line-feed” character (LF, ASCII code10), or the “EOI” (End Of Instruction) GPIB control line, or both.
- Instructions syntax is not case sensitive.

Instructions are divided into two groups:

1. **Commands** are used to configure the TUNICS-**Plus** by selecting operating modes and setting operating parameters such as laser output power, and more. (For example, to set the laser output power to 5 mW, user simply needs to type in the command string “P=5”.)
2. **Queries** are used to monitor parameter values at any time during system operation. (For example, the command string “P?” probes the level of laser output power. The TUNICS-**Plus** is addressed as the “talker” and returns the present parameter value, for example, “P=5” if the power output is 5 mW at the time of query.)

Each command/query is immediately handled by the TUNICS-**Plus** microprocessor. “Instruction_string” consists of a keyword with or without a question mark “?”. Some commands require a parameter. The parameter and the keyword are separated by a “=” sign or a “white space”.

Response messages

Upon completion of an instruction, TUNICS-**Plus** does not return any response message unless the command string is a query.

Value format

The following format rules apply to parameter values:

- Messages may contain one or more instructions. When a string containing several messages is sent, consecutive instructions must be separated by the semi-column “;” character.

- Numeric, non-integer values such as the laser output power (P) must be transmitted in a fixed-decimal-dot format (for example, P=5.1). However, integer values such as the validation byte of the request condition for service of the SRE instruction do not take a decimal-dot.
- White space¹³ characters are allowed in some places inside the strings sent to the TUNICS-**Plus**. They can be placed at the beginning of an instruction, after an instruction, as well as before, after, or in place of the “=” character. They cannot be inserted either in an instruction mnemonic and a numeric value, or between a mnemonic and the “?” character.
- Leading zero “0” characters are allowed at the beginning of a value.
- Digits after the decimal point may be sent or omitted.
- Decimal coma is not equivalent to decimal point.
- Values are ended with the “line-feed” character.
- Do not specify any units to parameter values. All values take the default units, factory-set for the TUNICS as shown in Table 6 - 11 below.

Table 6 - 11: Sample command syntax (GPIB)

Entry	Correct/Incorrect	Description
L= 1 530.2	Incorrect (no space allowed between "1" and "5")	Command error
L=1530.2	Correct	Set wavelength to 1530.200 nm
L=1530,2	Incorrect (no decimal coma allowed)	Command error
P=1 mW	Incorrect (unit abbreviations "mW" are not allowed)	Command error
P=1	Correct	Set power to 1 mW.
P=01	Correct	Set power to 1 mW.

Synchronization of Messages

To accelerate and secure the exchange of information between the computer/terminal and the TUNICS-**Plus**, we recommend using the value of the status word obtained by serial polling (see *Definition of the Status Word*).

This status word contains the information letting the computer/terminal know the status of commands received by the TUNICS-**Plus**, and therefore decide when a new instruction can be sent or when the requested parameter can be read.

Various flags in the status byte register (STB) indicate the outcome of the command last processed and the type of error event, if any (see also Table 5-8).

¹³

On the IEEE-488 interface, all characters with ASCII codes lower than or equal to 32, except the line-feed character (LF, ASCII code 10), are considered as white space characters.

Message synchronization is essential in the following two cases:

- **To ensure that a pending command has been completely executed before proceeding to perform other operations.** This can be checked through the OPC flag (bit 0 of STB). This flag remains false (0) as long as a command line has not been completely executed, and is set to true (1), afterwards. The computer should test this flag until it becomes true, and only then, proceed with the next instruction. In this case the STB byte must be read through a serial poll.
- **When a query statement has been sent, the computer must wait until the response is actually available before reading it.** This is signaled by the MAV flag (bit 4 of the STB byte). In this case the STB byte must be read through a serial poll. This flag remains false (its bit value is 0) until a complete message is available for reading. Therefore, the computer should constantly check this flag until it becomes true (bit value equal to 1), and only then, the response message associated with the query becomes available.

Error Handling

Two types of errors can occur in the transmission of instructions: command errors and value errors.

The system detects command errors each time a command string received by TUNICS-Plus is invalid due to a syntax error or to the use of an unknown mnemonic. The ERRC indicator in the status word is then activated.

A value error occurs when the command is valid, but the value of the parameter is incorrect, either because it could not be read or because it is outside the valid range. In both cases, the provided command is ignored and the ERRV indicator in the status word is activated.

Both ERRC and ERRV indicators are automatically reset when a new valid command is received.

Changing the GPIB Address

The TUNICS-Plus system GPIB link primary address is factory-set to 10. However, user can change the GPIB link address at any time by pressing the **2nd** + **Remote** keys. The current GPIB address is then displayed on the TUNICS-Plus screen, and a new value can then be entered. A valid GPIB primary address must be in the range of 1 to 30.

GPIB address change can also be performed in IEEE remote mode by typing the command: "GPAD=nn", where "nn" is an integer within the 1-30 range.

The new GPIB address is then stored in flash memory, and becomes the new default system setting. TUNICS-Plus retains the new GBIP address even after system turn-off.

Local Lockout

When the TUNICS-**Plus** receives a message via its IEEE-488.1 interface, it switches automatically to the remote operation mode, and the Remote LED on the front panel is lit. In remote operating mode, all other keys of the TUNICS front panel become inactive. TUNICS can return into local (or manual) operating mode by pressing the **Remote** key.

If the “Local lockout” instruction is provided via the GPIB interface, the local (or manual) mode can no longer be restored via the TUNICS front panel Remote key (pressing the Remote key will trigger a “Local lockout” message on the TUNICS display). To return into local mode, user must enter the GPIB command “**Go to local**”.

IEEE-488.1 Control Commands

IEEE Standard Commands

The following standard commands of the IEEE-488.1 interface can be processed by TUNICS-**Plus**:

Note: These correspond to basic commands whose syntax varies depending on the particular GPIB driver software you are using. For the exact syntax of standard commands, please refer to the user's manual that comes with your GPIB hardware and software.

Table 6 - 12: IEEE Standard commands

Function	Description
Remote change	Switches from remote to local mode or local to remote mode
Local lockout	Disables all controls from the front panel
Status reading	Conduct a serial poll.
Write	Write data to a device
Read	Read data from a device
Clear	Reset a specific device

Changing Parameter Settings

Table 6 - 13 lists the commands for each parameter setting of TUNICS-**Plus**. Each “n” represents a numeric value between 0 and 9.

Table 6 - 13: Parameter commands (GPIB)

Parameter Command	Meaning
I=nnn.n	Sets the laser current level in mA.
P=nn.nn ¹⁴	Sets laser output power to nn.nn mW.
P=±nn.nn ¹⁰	Sets laser output power to ±nn.nn dBm.
L=nnnn.nnn	Sets emission wavelength to nnnn.nnn nanometers.
F=nnnnnn.n	Sets emission optical frequency to nnnnnn.n GHz.
FSCL=nn.n	Switches to the Fine Scanning mode. (canceled by any other command). nn.n is the wavelength change in picometer (pm).

¹⁴

The units (dBm or mW) should be first selected with the DBM or MW command.

Table 6 - 13: Parameter commands (GPIB) (cont'd)

FSCF=n.nn	Switches to the Fine Scanning mode. (canceled by any other command). n.nn is the optical frequency change in GHz.
MOTOR_SPEED=nnn	Sets the sweep speed (in nm/s) used in swept-mode, hop-free scanning operation. Valid speed settings can range from 1 to 100 nm/s. Operational sweep speeds are: [1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 17, 18, 20, 22, 25, 29, 33, 40, 50, 67, 100] nm/s. Therefore, the speed setting you provided is automatically rounded to the nearest operational sweep speed. For example, if you provide "MOTOR_SPEED=060", the actual sweep speed implemented by the system is 67 nm/s. To check the current operational sweep speed, use the "MOTOR_SPEED?" query.
*SRE=nnn	Defines the conditions under which TUNICS-Plus will automatically send a service request (SRQ) to the computer/terminal.
GPAD=nn	Sets the GPIB address of TUNICS-Plus.
B_SUPPR=n	0=Eliminates the effect of the tuning micrometer screw backlash. 1=Recovers the backlash suppression procedure.
L_FEEDBACK=1	Enables the "active wavelength monitoring" while using an external wavemeter (see also page 46).
L_FEEDBACK=0	Disables the "active wavelength monitoring".
PCAL1 =nn.nn	Sets the calibration power (@ first wavelength).
PCAL2 =nn.nn	Sets the calibration power (@ second wavelength).

Changing Operating Mode

Table 6 - 14 lists the TUNICS-Plus operating mode commands.

Table 6 - 14: Operating mode commands (GPIB)

Command	Action
APCON	Turns on the Automatic Power Control.
APCOFF	Turns off the Automatic Power Control. Note: Using the APCON and APCOFF commands is the same as pressing the APC key in manual operating mode to toggle between constant power and constant current mode.
ACTCTRLON	Enables active cavity control.
ACTCTRLOFF	Disables active cavity control.
AUTO_CAL	Launches the system auto-calibration sequence. (Same as pressing the 2nd + I keys in manual mode.) Note 1: Using the AUTO_CAL command forces the system into performing an auto-calibration. Note 2: If the auto-calibration has failed, the message "referencing error" is displayed for a few seconds, and the TUNICS resumes normal operating status. Unlike the 2nd + I command that stalls the system when a referencing error occurs, there is no need to press Enter to restore the system.

Table 6 - 14: Parameter commands (GPIB) (cont'd)

DBM	Switches power units to dBm.
MW	Switches power units to mW.
ENABLE	Enables the laser output.
DISABLE	Disables the laser output.
CTRLON	Sets the linewidth coherence control ON.
CTRLOFF	Sets the linewidth coherence control OFF.

Note: Wavelength scans can only be performed manually, via the command keys at the TUNICS-Plus front panel.

System and Parameter Queries

Table 6 - 15 lists TUNICS-Plus system parameter queries. User can send queries to monitor system parameters in real time as well as prompt for current software version.

Note: A valid parameter value query consists of the command name immediately followed by the “?” mark.

Table 6 - 15: System parameter queries (GPIB)

Queries	Result
I?	Returns the value of the laser-diode current (mA).
P?	Returns the value of the optical output power (mW or dBm) ¹⁵ .
L?	Returns the value of the emission wavelength (nm).
F?	Returns the value of the optical frequency (GHz).
*IDN?	TUNICS-Plus returns the message “Anritsu,TUNICS_Plus,0,n.nn” where n.nn is the software release.
B_SUPPR?	Returns: "0" when backlash suppression is On; "1" when backlash suppression is Off.
MOTOR_SPEED?	Returns the operational sweeping speed. “MOTOR_SPEED=nnn” (in nm/s) ($1 \leq nnn \leq 100$)
ENABLE?	Returns the state of the laser output (ENABLE/DISABLE).
L_FEEDBACK?	Indicates whether the “active wavelength monitoring” function is in use. “1”: in use “0”: turned off “ERROR”: possible malfunction (See <i>Appendix A: Troubleshooting</i> for details.)
LIMIT?	Returns the status of the current limitation feature. “Yes” means the current is limited or “No” otherwise.
PCAL1?	Returns the calibration power at first wavelength LCAL1.
PCAL2?	Returns the calibration power at second wavelength LCAL2.

15

The units (dBm or mW) should be first selected with the DBM or MW command.
(See *Table 5-9: Operating mode commands*.)

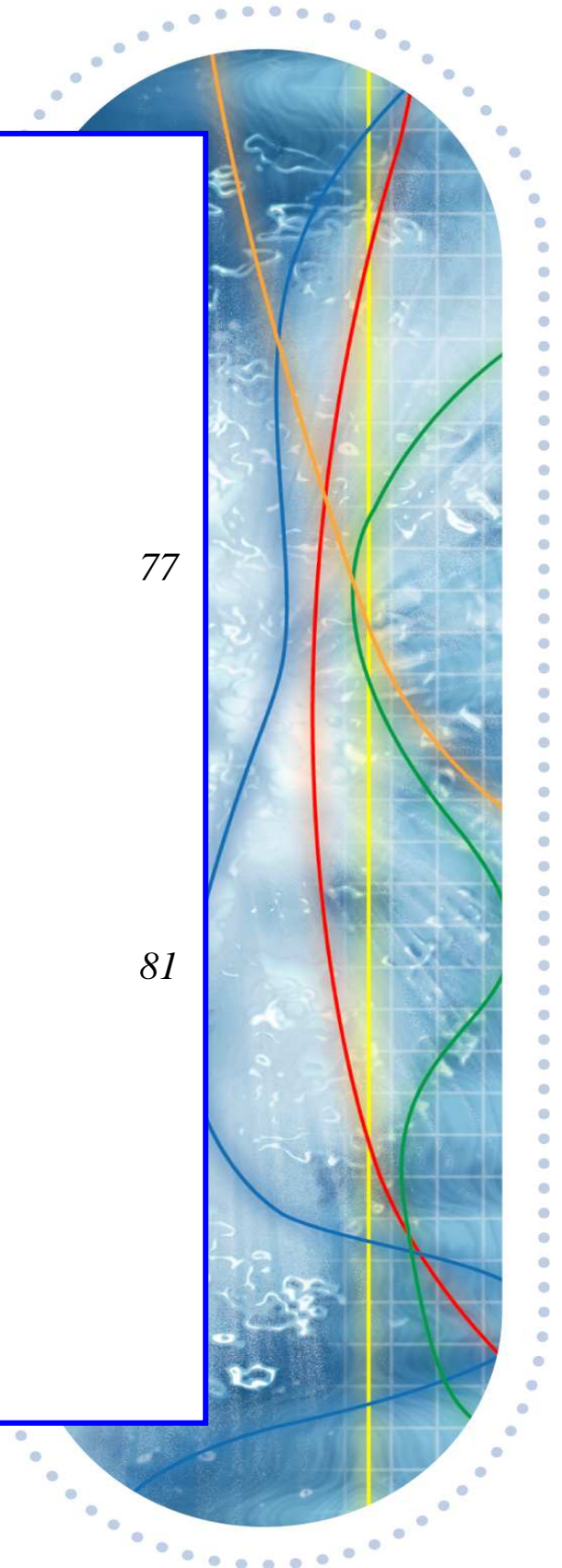
SERVICE INFORMATION

Performance Verification

77

Maintenance

81




Service Information

Performance Verification

Use the procedures in this section to verify the following warranted specifications of the TUNICS-**Plus** Tunable External Cavity Laser:

- Absolute wavelength accuracy
- Wavelength range and output power
- Tuning repeatability
- Side mode suppression
- Return loss

The recommended calibration interval is one year.

	WARNING	<p>To avoid eye damage, do not look into the laser output or into the end of an optical cable. Even though the laser radiation is not visible, the intense infrared light can cause eye damage.</p> <p>Always disable the laser output while connecting and disconnecting the fiber-optic cable.</p>
---	----------------	--

Equipment Required

Table 7 - 1 lists the equipment to carry out the performance verification procedures. The type and amount of connectors needed may vary depending on the particular equipment you use as well as the distance between the test equipment and the TUNICS-**Plus** tunable laser under test.

Table 7 - 1: Test equipment

Description	Minimum Requirements	Example Products
Optical power meter	Wavelength range: 1260 nm to 1640 nm	RIFOCS 578L ANDO AQ 2140
Optical Spectrum Analyzer (OSA)	Wavelength range: 1250 nm to 1650 nm	Yenista OSA425
Wavemeter	Wavelength range: 1260 nm to 1650 nm	BURLEIGH WA1000
Jumper	FC/APC to FC/PC	
Inline optical adapter	FC to FC	

Preparation

Warm up the TUNICS-**Plus** tunable laser under test and the test equipment for 2 hours at an ambient temperature (between 20° and 30° C).

Initializing the System

To ensure safe operation, use the following procedure to initialize the system:

- Make sure the key on the front panel is in the Standby position.
- Connect the power supply cable to the rear panel of the TUNICS-**Plus** unit on one end and to the proper power outlet on the other end.
- Set the rear panel MAINS switch to ON.
- Turn the front panel on/off key to the right. During the initialization phase, the message “**Initializing...**” is displayed.
- The message “**Referencing...**” is then displayed to indicate that the system performs its self-calibration (internal referencing of the wavelength). This message clears once the TUNICS-**Plus** tunable laser is initialized and ready for operation.
- At the end of the initialization routine, the tunable laser is ready to enable. The message “**disabled...**” appears on the display.

	IMPORTANT	Once the initialization process is complete, the laser is still not enabled for safety reasons, and the message “disabled...” is displayed. To enable the tunable laser, user must press the Enabled key on the TUNICS- Plus front panel.
	CAUTION	Avoid switching the unit off before the system is fully initialized.

Connecting the Laser Output

To connect the laser output, remove the protective cap and use an optical fiber with a FC/APC type connector.

	CAUTION	To avoid damage to the laser output, always use a FC/APC type connector. To connect other types of connectors such as FC/PC, make sure to use the appropriate FC/APC converter. To help protect the laser output, leave a fiber-optic cable connected to the unit while making connections to other devices. When the laser output is not connected, always install the protective cap.
--	----------------	--

To achieve optimum system performance, optical connections must be clean and free of contamination. Cleaning instructions for optical connections can be found in the *Maintenance* section in this chapter.



Enabling the Laser Output

For safety reasons, the laser output is enabled only when user presses the **Enabled** key.

- To enable the laser output, press the Enable key. The Enable LED is lit to indicate that the TUNICS-**Plus** is emitting a laser output.
- To disable the laser output, simply press the Enable key again and check that the Enable LED went off.

Shutting Down the Laser

To shut down the laser, turn the front panel key counterclockwise to the Standby position. The message “**Parking...**” is displayed and the optical head is moved automatically to its parking position. Wait a few seconds till the system shuts down by itself.

	WARNING	Do not unplug the fiber-optic cable from the laser output connector while the laser is still enabled, as exposure to the laser beam is potentially harmful.
	CAUTION	Follow these guidelines while shutting down the system: <ul style="list-style-type: none"> ▪ Do not turn on and off the system repeatedly and too rapidly, and particularly after turning off the laser, wait at least thirty seconds before initializing the laser again. ▪ Avoid switching off the unit before the system is fully initialized.

Absolute Wavelength Accuracy

This procedure verifies that the wavelength settings agree with the measured wavelengths over the entire specification range. If these values do not match the absolute wavelength specification, user may calibrate the lower and upper wavelength setting limits by using an external wavemeter, as described below:

- Connect the laser to the wavemeter.
- Set the output power to 0 dBm using the front-panel power controls: (Press P > 0 > ENTER).
- Set the wavelength of the tunable laser to the lower limit of the wavelength range according to the **TUNICS-Plus** model (see *Appendix B: Valid Range for Parameter Settings* in the *Appendices* chapter).
- Enable the laser output and verify the wavelength at the lower limit ± 0.04 nm.
- Repeat the check for each increment of 5 nm from the lower limit to the upper limit listed in *Appendix B*.
- Verify that the measured wavelength is no greater than ± 0.04 nm away from each wavelength setting.

Wavelength Range and Output Power

This procedure verifies that the laser output can achieve the rated power at the specified limits of the wavelength range. If the power settings do not exactly match the measured power, you may calibrate the power settings from the front panel, as follows:

- Disable the laser output.
- Connect the laser to the optical power meter.
- Set the wavelength of the tunable laser and the power meter to the lower limit of the wavelength range listed in Table 2 - 1 in the *General Information* chapter.
- Enable the laser output and adjust the power setting as necessary to verify that the laser can achieve the rated output power at the lower limit of the wavelength range.

- Set the wavelength of the tunable laser and the power meter to the upper limit of the wavelength range listed in Table 2 - 1.
- Enable the laser output and adjust the power setting as necessary to verify that the laser can achieve the rated output power at the upper limit of the wavelength range.

Side Mode Suppression

This procedure verifies the suppression of any wave emissions that are more than 1 nm away from the selected wavelength. This procedure requires the use of an optical spectrum analyzer:


- Connect the laser output to the optical spectrum analyzer.
- Set the laser wavelength at mid-range of the wavelength specification range listed in Table 2 - 1 in the *General Information* chapter.
- Set the laser output power to 0 dBm
- Configure the optical spectrum analyzer as follows:

Control	Setting
Center wavelength	Same as laser wavelength
Span	20 nm
Resolution	0.1 nm
Reference level	0 dBm
Level scale	10 dB/div
Scan	Repeat

User can now verify that at wavelength deviations larger than 1 nm from the centered wavelength of the laser peak, the power level of any residual emissions should be lower than -45 dB down from the peak level.

Wavelength and Power Calibration

During the course of the performance verification it is possible that the settings for wavelength and power may not agree with the measured values. If this occurs, you may calibrate wavelength and power settings from the front panel.


	CAUTION	The procedures to calibrate the wavelength and power settings erase factory-set calibration data. Do not perform these procedures unless you have calibrated sources.
---	----------------	---

To perform the wavelength and power calibration procedures, refer to the *Wavelength Calibration and Power Calibration* sections in the *Operation* chapter.

Maintenance

User maintenance of the TUNICS-**Plus** is limited to the following:


- Cleaning optical connectors for optimum power
- Checking the performance of optical connections and cables
- User calibration
- Replacing the line fuse

	WARNING	<p>DO NOT MAKE ANY SERVICE OR MAINTENANCE TO THE TUNICS-Plus LASER HEAD:</p> <p>Under no circumstance must the user make any service or maintenance repairs of any kind to the TUNICS-Plus laser head. The user must not remove the protective cover of the TUNICS-Plus laser head. Refer servicing to authorized Yenista personnel only.</p>
---	----------------	---

Cleaning and Caring for Optical Connectors and Fibers

- Handle optical fiber with appropriate care and keep optical connectors free of outside contaminant to preserve optical-signal integrity.
- Use the following items to clean the optical connectors:
 - clean compressed air
 - fiber-optic cleaning swabs
 - isopropyl alcohol
- Follow these steps to clean the optical connectors:
 1. Hold the can of compressed air upright and spray the can into the air to purge any propellant.
 2. Spray the clean compressed air on the connectors to remove any loose particles or moisture.
 3. Moisten a clean optical swab with isopropyl alcohol, then lightly swab the surfaces of the connectors.
 4. Spray the clean compressed air on the connectors again to remove any loose particles or isopropyl alcohol.

Note: Cleaning kits for optical connectors are available from many fiber optic suppliers.

	CAUTION	<p>To optimize the performance of the system and to prevent loss of optical power or damage to the optical connectors, keep the connectors clean at all times.</p> <p>When cleaning the connectors with a swab, use gentle circular motions. Use only high quality cleaning supplies that are non-abrasive and leave no residue.</p> <p>To reduce the need for cleaning, immediately replace protective caps on the optical connectors when not in use.</p> <p>Do not over stress or sharply bend optical fiber beyond tolerances specified by the manufacturer.</p>
---	----------------	--

- To clean the internal connector, unscrew the plug, slightly pull it towards you, and clean the internal connector. Then put every thing back in place. This operation must be made carefully in order not to damage the internal connector near the fiber cable.

Cleaning of the Instrument

If necessary, follow these steps to clean the TUNICS-**Plus** unit:

- Unplug the power cord from the mains supply before cleaning.
- Clean the case with a slightly damp cloth with an isopropyl alcohol liquid.

	CAUTION	Do not use chemically active or abrasive materials to clean the instrument.
--	----------------	---

User Calibration

To externally verify the performance of the TUNICS-**Plus**, refer to the procedure described in the *Wavelength Calibration and Power Calibration* sections in the *Operation* chapter. If the performance does not meet warranted specifications, a complete calibration of the instrument may be required.

Some calibration steps require the adjustment of internal components including the laser head. Therefore, a complete calibration must be performed by authorized Yenista personnel only.

Replacing the Line Fuse

Before replacing the line fuse, be sure to unplug the power cord.

- Insert a small flat-blade screwdriver into the notch just inside the power cord socket of the line filter (see Figure 7 - 1).
- Use the notch to pull the fuse holder straight out to remove the fuse.
- Replace the fuse (4.0A 250V SLOW/(T)) in the fuse holder and snap the fuse holder back into the lined filter.

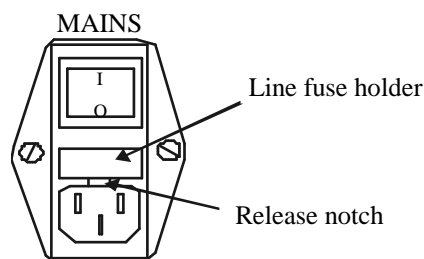


Figure 7 - 1: Replacing the line fuse

Packaging for Shipment

When you receive the instrument, please keep the original packaging for potential use.

Use the original packaging in case you need to return the TUNICS-**Plus** tunable laser source unit to Yenista for servicing or calibration. For instructions on returning the instrument, please contact Yenista (see *Contact Information* on page 8).

APPENDICES

<i>Appendix A:</i> <i>Error and Warning Messages</i>	85
<i>Appendix B:</i> <i>Valid Range for Parameter Settings</i>	86
<i>Appendix C:</i> <i>Certification and Compliance</i>	87
<i>Appendix D:</i> <i>Replaceable Parts</i>	88
<i>Index</i>	89

Appendices

Appendix A: Error and Warning Messages

Message	Possible cause	Troubleshooting
WARNING: User calib. Err press “enter”	User calibration is not correct.	As the user calibration is not valid, TUNICS- Plus implements the default calibration value instead.
ERROR: No factory calib. Turn Power off	No factory calibration.	Turn off the system and check the integrity of electrical connections. If this error message shows continuously, please contact a service engineer immediately.
VALUE error Press Enter.	The system parameter you are trying to configure is outside the valid range. (See <i>Appendix B: Valid Range for Parameter Setting</i>)	Press Enter. All changes are discarded and the previous parameter value is restored.
Limit Switch Turn Power off	Mechanical problem.	Turn off the system and contact a service engineer immediately.
EEPROM error Turn Power off	Electrical problem.	Turn off the system and contact a service engineer immediately.
Referencing error. Press Enter	The system's internal wavelength referencing operation has failed.	Press Enter to restore the current system settings. Internal wavelength referencing remains unchanged.
λ FEEDBACK:ERR	Active wavelength monitoring does not function properly. Possible causes are: 1. No communication between the TUNICS- Plus tunable laser and the wavemeter. 2. The wavemeter may not be compatible with the SCPI command set. 3. The wavemeter remote interface may not be configured properly. 4. No optical signal. 5. Wavemeter malfunction.	Press Enter to restore the current system settings. Internal wavelength referencing remains unchanged. 1. Check the RS-232 serial connection between the two devices. 2. Please contact a Yenista representative for a list of compatible wavemeters. 3. Configure the wavemeter serial port with the following parameters: <ul style="list-style-type: none"> • RS-232 → ON • Baud rate → 9600 • Handshaking → None 4. Press the TUNICS- Plus front-panel Enable key to turn on the laser output (the “Enable” LED is lit).

Appendix B: Valid Range for Parameter Settings

TUNICS-Plus Model	Plus O	Plus E	Plus S
Laser diode current level (I)	See attached test report	See attached test report	See attached test report
Beam output power (P)	0.2 – 10 mW	0.2 – 10 mW	0.2 – 10 mW
Upper and lower wavelength scan limits (nm)	1240-1380	1330-1470	1410-1550
Wavelength scan steps (nm) <ul style="list-style-type: none"> • Step-mode • Swept-mode 	0.001-150 continuous	0.001-150 continuous	0.001-150 continuous
Time at each scan step ¹⁶ (Pause Timer)	0.1 - 25 sec	0.1 - 25 sec	0.1 - 25 sec

TUNICS-Plus Model	Plus S WB	Plus CL	Plus CL WB
Laser diode current level (I)	See attached test report	See attached test report	See attached test report
Beam output power (P)	0.2 – 10 mW	0.2 – 10 mW	0.2 – 10 mW
Upper and lower wavelength scan limits (nm)	1370-1560	1500-1640	1470-1660
Wavelength scan steps (nm) <ul style="list-style-type: none"> • Step-mode • Swept-mode 	0.001-150 continuous	0.001-150 continuous	0.001-150 continuous
Time at each scan step ¹⁶ (Pause Timer)	0.1 - 25 sec	0.1 - 25 sec	0.1 - 25 sec

¹⁶ Step-mode scanning only.

Appendix C: Certification and Compliance

Category	Standards or Description
EC Declaration of Conformity-EMC	Compliance was demonstrated to the following specification as listed in the Official Journal of the European Union: Directive 89/336/EEC for Electromagnetic Compatibility. EN 50081-1 Emission: EN 55022 Class B Radiated and Conducted Emissions EN 50082-1 Immunity: EN 61000-4-2 Electrostatic Discharge Immunity IEC 60801-3 RF Electromagnetic Field Immunity EN 61000-4-4 Electrical Fast Transient/Burst Immunity
EC Declaration of Conformity-Low Voltage	Compliance was demonstrated to the following specification as listed in the Official Journal of the European Union: Low Voltage Directive 73/23/EEC, amended by 93/68/EEC. EN 61010-1: 1993/A2: 1995 Safety requirements for electrical equipment for measurement control and laboratory use. EN 60825- 1/A1: 1994 Safety for Laser Products, Part 1. Equipment classification requirements and user's guide.
U.S. Nationally Recognized Testing Laboratory Listing	UL3111-1 Standard for electrical measuring and test equipment.
Canadian Certification	CAN/CSA 1010-1-92 including MOD: 1997 Safety requirements for electrical equipment for measurement, control, and laboratory use.

Safety Certification and Compliance

Equipment Type	Test and Measuring
Safety Class	Class 1, as defined in IEC-61010-1 annex H-Grounded product.
Mains supply voltage	100-240 V AC with $\pm 10\%$ fluctuation max.
Installation (overvoltage) category	Overvoltage category II, as defined in IEC 61010-1, annex J. Local-level mains (wall sockets). Equipment at this level includes appliances, portable tools, and similar products. Equipment is usually cord-connected.
Pollution Degree	Pollution Degree 2, as defined in IEC 61010-1. Note: rated for indoor use only. Normally only dry, nonconductive pollution occurs. Occasional and temporary conductivity caused by condensation may occur. This location is a typical office/home environment. Temporary condensation occurs only when the product is left unused.
Safe operating temperature range	+5°C to +40°C.
Maximum relative humidity	80% for temperatures up to 31°C decreasing linearly to 50% relative humidity at 40°C.
Altitude (maximum operating)	2000 meters.

Appendix D: Replaceable Parts

For information about replaceable parts, contact your Yenista sales representative.

Standard Accessories

Name	Description
Fuse	4A T 250V
Manual, tech: instruction	TUNICS- Plus instruction manual
PM optical fiber cable	Single-mode polarization maintaining cable
Optical fiber cable	Single-mode cable
Power key	Front panel keylock
Rack-mount handles	2 rack-mount handles

Optional Power Cords

Name	Description
North America 110 V, 60 Hz	CSA STD C22.2 NO21-M1984, ECN 436C approved, 5-15P attachment plug cap molded onto No18 AWG, type SVT,SJT flexible cord, 2 meter (78.7 inch), 125V/10A AC, US
Universal European 220 V, 50 Hz	DIN 49441/2, CEE7, DIN 0625, VDE approved, IEC attachment plug cap, HO5VV-F3G 0.7 mm ² flexible cord, 2.5 meter (98.4 inch), 250V/10A AC, EUROPEAN
United Kingdom 240 V, 50 Hz	BS4491.EN60.320,BASEC BS6500, ASTA BS 1363A, VDE, SEV, BSI approved, 1.0 mm ² flexible cord, 2 meter (78.7 inch), 250V/10A AC, 25.4 mm fused UK plug (13A fuse)

Index

- Absolute Wavelength Accuracy, 79
- Accessories, 20
 - Replaceable Parts, 4, 88
 - Standard Accessories, 4, 19, 88
- Active Cavity Control, 13, 45
 - Disable/Enable, 28, 45
 - Operating requirements, 45
- Active Wavelength Monitoring, 46
- APC function key, 37
- Appendices, 4, 85
- Auto-calibration, 28
- Calibration
 - Calibration interval, 77
 - Certificate of traceable calibration, 19
 - Power calibration, 40
- Certification and Compliance, 87
- Coherence Control, 29, 42, 55
- Contact Information
 - Service Information, 4, 77
- Current
 - Diode current, 37
- Current limitation, 67
- DATA Area, 27
- Diode Current
 - Cutoff current, 51
 - Saturation current, 51, 54
- Error and warning messages, 4, 85
- Front panel, 27
- FSC See *Fine Scanning Mode*, 28, 29, 31, 37, 50, 55
- Functions (*Advanced*)
 - Backlash Suppression, 42
 - Coherence Control, 29, 42
 - GPIO Address, 42, 71
 - Power Calibration, 40, 80, 82
 - Step-Mode Wavelength Scanning, 38
 - Wavelength Calibration, 80, 82
 - Wavelength Sweep
 - Sweep speed, 40
- Functions (*Most Commonly Used*)
 - Diode Current, 37
 - Emission Wavelength, 36
 - Fine Scanning Mode, 37
 - Optical Output Power, 36, 38
 - Switching between constant-power and constant-current modes, 37
 - Wavelength Step, 36
- GPIO
 - GPIO Address, 71
 - GPIO Link, 66
 - GPIO Standard, 66
- High Frequency Modulation, 52, 53
- IEEE-488.1
 - Connector, 31
 - Remote control mode, 66
- INIT, 63
- Initializing, 22, 78
- Inputs/Outputs
 - λ , I, P, Sync (Monitor outputs), 31, 55
 - FSC (Wavelength modulation input), 55
 - IEEE-488.1 Interface, 50
 - Mod Out (Internal modulation output), 55
 - Mod Sync. (Synchronization signal output), 55
 - Remote Interlock, 6, 31, 56
 - RS-232 C Interface, 50
 - Sync. (Motor Sync.), 31, 50, 55
- Installation, 21
- Internal Wavelength Referencing, 13, 43
 - Auto-calibration, 28, 43
 - Self-calibration at start-up, 43
- Laser
 - Absolute maximum ratings, 6
 - Class, 6
 - Enable, 27, 63, 74
 - Initializing, 22
 - Optical connection, 22
 - Shutting down, 23, 79
 - Specifications, 6
- Lim, 36, 39
- Linewidth See Specifications, *See also*
 - Coherence Control, 14, 42, 55, 64, 74
- Local Lockout See *Operating mode*, 72
- Low Frequency Modulation, 50, 51
 - I_{diode} , I_{bias} , I_{in} , I_{cutoff} , I_{max} , 50, 51
- Maintenance, 81
- Message
 - Errors and Warnings, 4, 85
 - Message available (MAV), 67, 68
 - Message Format (RS-232 C), 60
 - Response Messages (RS-232 C), 61
 - Synchronization (GPIO), 70
 - Syntax and Message Conventions (GPIO), 69
- MODE Area, 28
- Mode-locked operation, 14, 53
- Models
 - Ordering codes, 19
- Modulation
 - High Frequency Modulation, 52
 - Internal Modulation, 55
 - Low Frequency Modulation, 50
 - Synchronization Signal, 55
 - Wavelength Modulation, 55

Appendices

- Operating mode
 - Local Lockout, 72
 - Manual mode, 3, 28, 35
 - Operating mode commands (GPIB), 73
 - Operating mode commands (RS-232 C), 63, 64
 - Remote control, 4, 59, 66
- Operating Mode
 - Constant-Current mode, 28, 37
 - Constant-Power mode, 28, 37
- Optical connection, 22
- Optical design, 12
- Options, 4, 19
 - Ordering codes, 20
- Ordering codes, 20
- Ordering information, 19
- Parameter
 - Command parameters (GPIB), 72, 73, 74
 - Command parameters (RS-232 C), 62
 - Reading system parameters (RS-232 C), 64
 - Setting (General Principle), 29
 - Setting (RS-232 C), 62
 - System parameter queries (GPIB), 74
 - Valid range*, 36, 86
- Parameter settings
 - Changing parameter settings, 32
 - Valid range*, 4, 86
- Parking, 23, 79
- Pause Timer, 38, 39, 86
- Power
 - Absolute maximum ratings, 6
 - Optical output power, 36
 - Standard, 6
- Query, *See Synchronization of messages*
 - System parameter queries (GPIB), 74
 - System parameter queries (RS-232 C), 64
- R option* *See Specifications*, 15
- Rear panel, 31
 - Connectors, 31, 50
- Remote control, 16
 - IEEE-488.1, 66
 - RS-232 C, 59
- Remote Interlock, 6, 31, 56
- Resolution
 - Wavelength setting resolution, 13
- RIN *See Specifications*, 15
- RS-232 C
 - Connector, 31
 - Remote control activation, 42
 - Remote control mode, 59
- Safety
 - Electrical Safety, 7
 - Equipment Safety Warning, 7
 - Laser Safety, 6
- Scan
 - Fine scanning mode, 28, 37
 - Pause Timer, 38, 39, 86
 - Scan interval, 38
 - Scan step, 38, 39, 86
 - Step-Mode Wavelength scanning, 38
 - Stop scan, 39, 40
 - Wavelength scan limits, 38, 39, 86
- Serial polling, 67, 68, 70
- SETUP Area, 29
- Side Mode Suppression, 15, 80
- Side Mode Suppression Ratio *See Specifications*, 15
- Specifications, 15
- Status Word
 - Definition of the Status Word, 67
 - Use of the Service Request (SRQ), 68
- Sweep *See Wavelength Sweep*
- System
 - Initialization, 22
 - Installation, 21
 - Shutting down, 23, 79
- Tuning repeatability *See Specifications*, 15, 77
- Units
 - dBm, mW, 27, 63
 - nm, GHz, 27
- Unpacking, 19
- Value error (ERRV), 67, 68, 71
- Wavelength accuracy *See Specifications*, 15, 77
- Wavelength range *See Specifications*, 6, 15, 77
- Wavelength setting resolution* *See Specifications*, 15
- Wavelength Sweep, 39
- wide-band option (/WB), 20
- Wide-band option (/WB) *See Specifications*

